

# Mobicols - Columns


Tools for diverse Lab Applications

SEPARATION

ISOLATION

TREATMENT

PURIFICATION

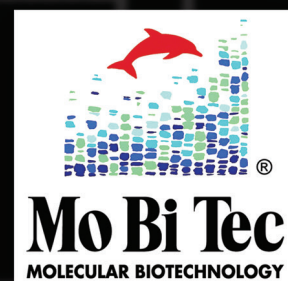





## Chapter I

**Mobicol – A Practical Tool for every Lab**

Page 3

## Chapter II

**Spin Columns for your own Matrices or  
other Compounds**

Page 11

## Chapter III

**Matrix-filled MobiSpin Columns for  
Purification of Nucleic Acids**

Page 15

## Chapter IV

**Mobicols for Affinity Purification of  
Biomolecules**

Fractionation and Purification of Antibodies  
Immobilized Metal Ion Affinity Chromatography

Page 21

## Chapter V

**Mobicols with Enzymatically Active Matrices**

Immobilized Proteinases  
Immobilized Nucleases  
Other Immobilized Enzymes

Page 25



## Introduction

Mobicols can be used to make your own spin/affinity columns with the matrix of your choice (**Chapter II**).

Used as a Mini Column they feature Luer-lock connections on top and bottom. Using the Luer-lock adapter, a syringe can be connected to the top of the Mobicol as a reservoir for sample or buffer application. To enable a connection to larger reservoirs, a Luer-adapter is available to attach tubing to the top cap. As Spin Columns, Mobicols are placed in a 1.5 ml or 2 ml microtube and can be centrifuged in a microcentrifuge.

Spin Columns can be used for the purification of biomolecules from small contaminants (e.g. desalting) (**Chapter III**).

Mobicols filled with a Protein A or Protein G matrix are useful for the purification of antibodies from cell culture supernatant, plasma or serum. Prefilled Mobicols are available in Starter Kits with all required buffers, columns and accessories (**Chapter IV**).

Applying Mini Columns, you can take advantage of highly active immobilized enzymes, packed into small Compact Reaction Columns (CRC) (**Chapter V**).







# Chapter I


## **Mobicol – A Practical Tool for every Lab**





## Mobicol "F"

Mobicols are plastic columns with particular handling advantages:

- Mobicols are compatible with laboratory standards
- Mobicols can be centrifuged in microcentrifuges
- Mobicols have a Luer-lock adapter
- Smallest and largest volumes can be treated easily

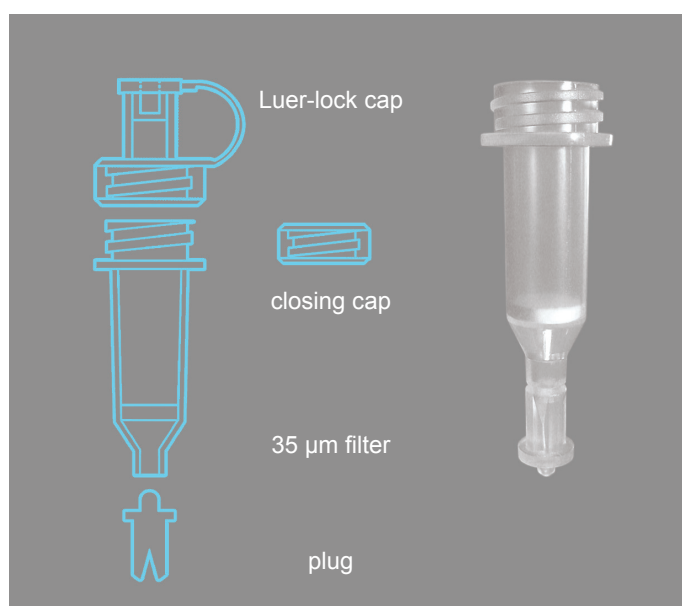

### Mobicol "F" with 35 µm filter in 2 ml collection tube

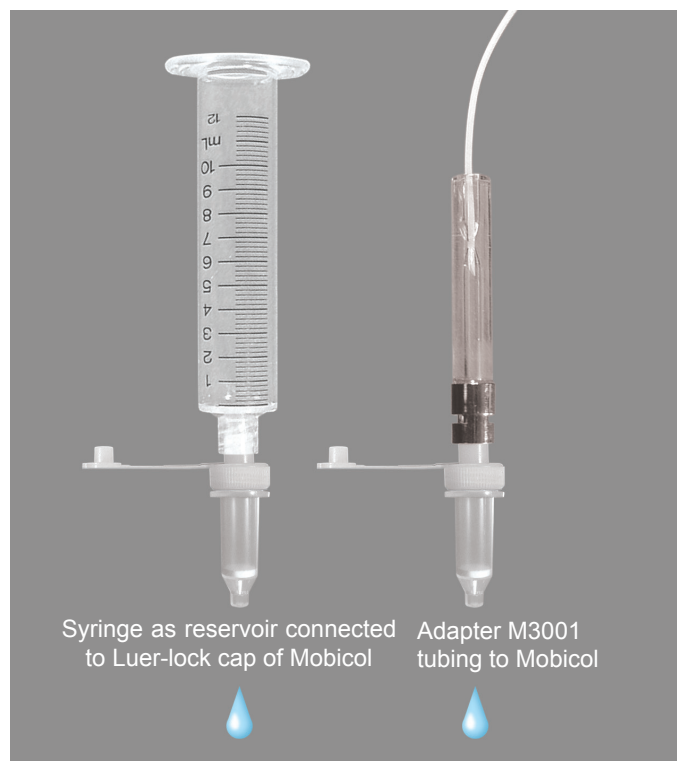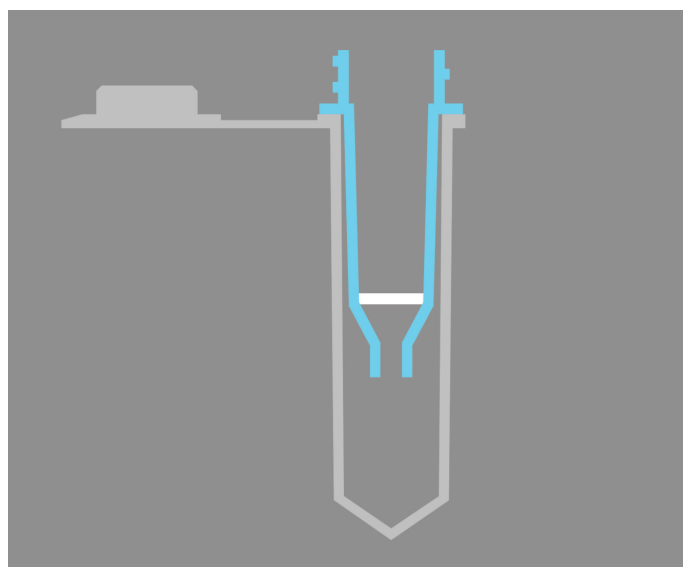



Filter of 35 µm pore size is inserted. A separately available upper filter (Order No. M523515) conveniently caps the resin bed.

Mobicol "F" can be closed tightly with a screw-on cap and a snap-off plug at the outlet.

Mobicol "F" can be autoclaved at 110 °C for 10 minutes.

### Mobicol "F" can also be used as Spin Column

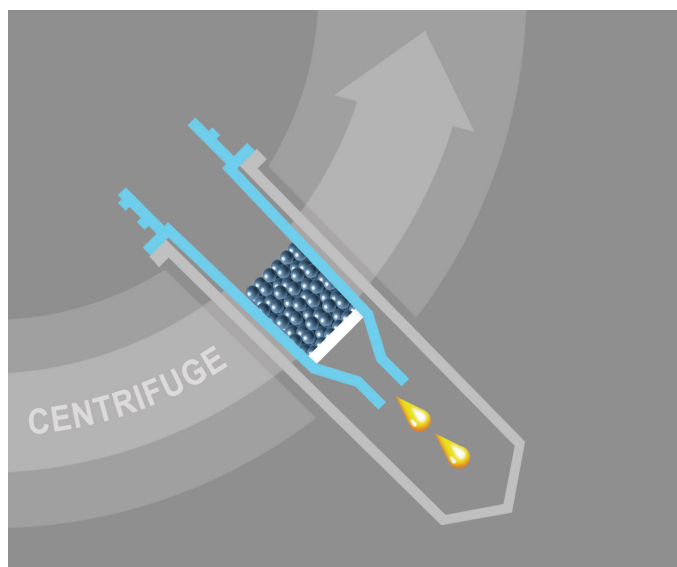

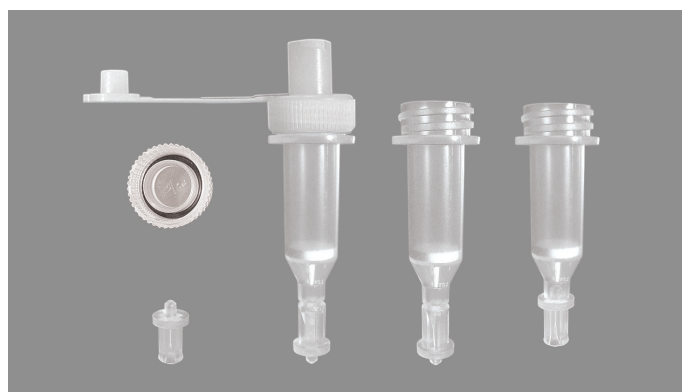

Mobicol "F" (M105035F) with fixed snap-off plug for the outlet is supplied with 2 caps: a Luer-lock cap and a screw cap. The Mobicol "F" can be closed tightly with screw cap and plug. A filter of 35 µm pore size is inserted; an additional upper filter (35 µm) is available.

## Applications

Use it with affinity matrix for:

- **Purification of tagged proteins**
- **Purification of peptides**
- **Purification of native proteins** using immobilized IgG, NTP-binding proteins (e.g. kinases), amino acid binding proteins, β-lactamase, Igase, peroxidase, EPO (erythropoetin), β-galactosidase
- **Fractionation and purification of antibodies** using protein A/G, immobilized epitopes, peptide matrix
- **Purification of glycoproteins and carbohydrates**
- **Purification of fatty acids and fatty acid binding proteins**
- **Whole mount embryo *in situ* hybridization**

## Mobicol "F" is approved

- Hundreds of satisfied customers are using Mobicols already for many years
- Many of referenced applications

## Mobicol "F" is versatile

- Can be used in flow-through mode or as Spin Column
- Useful for size exclusion or affinity purification of compounds like nucleic acids and proteins (native or tagged)
- Compatible with small and large sample volumes
- Upper filters of 35 µm pore size are separately available
- Screw cap with rubber seal
- Mobicol can be filled and stored or shipped without leakage
- Reverse side of snap-off plug can be used for resealing
- Empty columns for other filter pore sizes available

## Mobicol "F" is compatible with lab standards

- Can be centrifuged in standard microcentrifuge
- Luer adapters are available for inlet and outlet tubing
- Luer-lock cap is compatible to Luer syringes
- Autoclavable (with filters at 110 °C) for 10 min

## Order Information

Order No.	Description	Amount
<b>M105035F</b>	Mobicol "F" with fixed outlet plug, inserted 35 µm filter and 2 different screw caps	50
<b>M1050</b>	Mobicol "F" with fixed outlet plug and 2 different screw caps, without filters	50
<b>M1053</b>	Mobicol "F" with fixed outlet plug and screw cap, without filter	50
<b>M2210</b>	Filter (large) 10 µm pore size, for Mobicol M1002, M1003, M1050 & M1053	50
<b>M523515</b>	Filter (large) 35 µm pore size, for Mobicol M1002, M1003, M1050 & M1053	50
<b>M2290</b>	Filter (large) 90 µm pore size, for Mobicol M1002, M1003, M1050 & M1053	50


For column accessories see page 7.

Also available prefilled with immobilized enzymes: immobilized proteases, nucleases, alk. phosphatase, or β-galactosidase (Chapter V).



## Mobicol "Classic"

**Use for Affinity Chromatography (M1002)  
with inserted filters (large)**

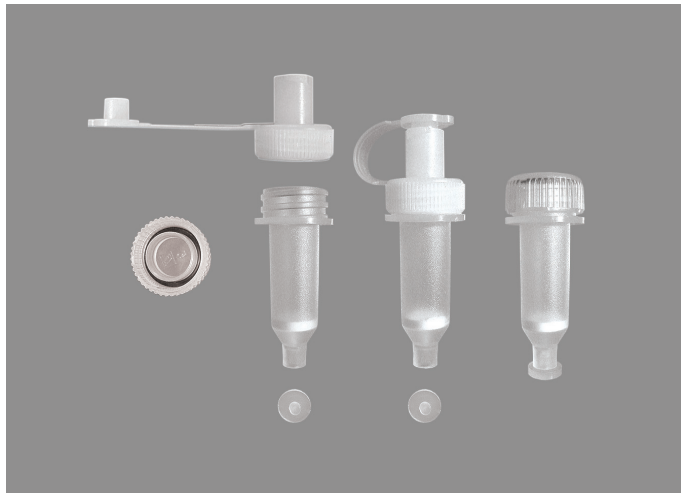

**Use as Spin Column (M1003)  
with inserted filters (small)**



The Mobicol "Classic" is supplied with one or two screw caps (Luer-lock cap and screw cap) and a plug for the outlet. Small and large filters with different diameters and pore size (10  $\mu\text{m}$ , 35  $\mu\text{m}$  and 90  $\mu\text{m}$ ) are separately available, so you can use the filter of your choice. With these features the Mobicol "Classic" can be used for a variety of applications.

**Ask for  
Free Samples!**

Filters and insertion tools



## Mobicol "Classic"

Order No.	Description	Amount
<b>M1002</b>	Mobicol "Classic" with 2 different screw caps, without filters	50
<b>M1003</b>	Mobicol "Classic" with 1 screw cap, without filters	50
<b>M2110</b>	Filter (small) 10 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
<b>M513515</b>	Filter (small) 35 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
<b>M2190</b>	Filter (small) 90 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
<b>M2210</b>	Filter (large) 10 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
<b>M523515</b>	Filter (large) 35 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
<b>M2290</b>	Filter (large) 90 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50



## Accessories for Affinity Chromatography Columns

The polypropylene Luer adapter (M) (M3003) for tubing connection to the inlet of the column allows a flow-through mode at low pressure, whereas the Luer adapter (F) (M3004) is suitable for a tubing connection to the outlet of the column. Furthermore, the Luer adapter (M/M) (M3007) facilitates the connection of the inlet to a Luer-lock system.

The Tubing-to-Luer-lock adapter (M3001) enables a safe connection from a pump to the inlet of a Mobicol or a 2.5 ml column.

The one-way stopcock, Luer-lock (F) (M3006) is a convenient tool for controlling flow-through at the in- or outlet.



**M3003**  
Luer adapter (M) for  
Luer-lock cap



**M3004**  
Luer adapter (F) for outlet



**M3001** Tubing-to-Luer-lock adapter (connected to Mobicol "Classic" (M1002) with inserted large filter; not enclosed)



**M3006**  
Stopcock, one-way,  
Luer-lock (F)



**M3007**  
Luer-lock adapter (M/M)


## Column Accessories

Order No.	Description	Amount
<b>M3001</b>	Tubing-to-Luer-lock adapter (material outside: metal; inside: glass and PTFE; with 1 m PTFE tubing) for Mobicols or laboratory columns	1
<b>M3002</b>	Luer adapter, (M+F) 10 each	20
<b>M3003</b>	Luer adapter for the top-cap (M)	20
<b>M3004</b>	Luer adapter for outlet (F)	20
<b>M3007</b>	Luer-lock adapter (M/M)	1
<b>M3005</b>	Bottom plugs for Mobicol M1002	50
<b>M3006</b>	Stopcocks, one-way, Luer-lock (F)	4
<b>M3009</b>	Luer-lock caps	50



## Related Products

For volumes larger than 1 ml (Mobicols), 2.5 ml, 5 ml and 10 ml columns are available with filters inserted. These syringe based columns have Luer-lock adapters at their in- and outlet, so that they have the same handling advantages as Mobicols.

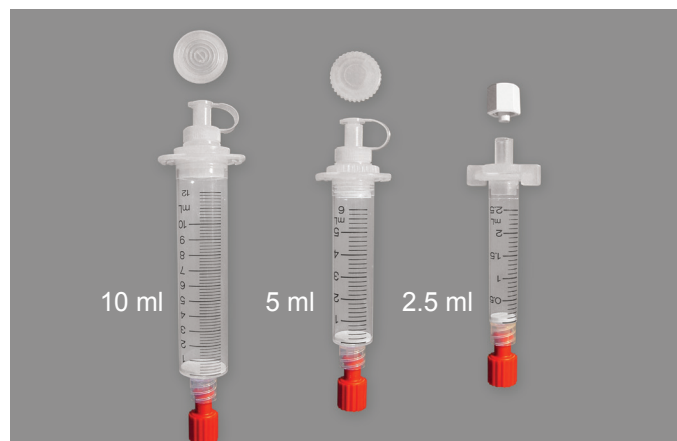

## Laboratory columns

Order No.	Description	Amount
<b>S10121</b>	2.5 ml columns (10 µm filter pore size)	20
<b>S1012</b>	2.5 ml columns (35 µm filter pore size)	20
<b>S10129</b>	2.5 ml columns (90 µm filter pore size)	20
<b>S10131</b>	5 ml columns (10 µm filter pore size)	20
<b>S1013</b>	5 ml columns (35 µm filter pore size)	20
<b>S10139</b>	5 ml columns (90 µm filter pore size)	20
<b>S10141</b>	10 ml columns (10 µm filter pore size)	20
<b>S1014</b>	10 ml columns (35 µm filter pore size)	20
<b>S10149</b>	10 ml columns (90 µm filter pore size)	20
<b>S1210</b>	Upper filter 10 µm for 2.5 ml columns	20
<b>S1235</b>	Upper filter 35 µm (thickness 1.5 mm) for 2.5 ml columns	20
<b>S123532</b>	Upper filter 35 µm (thickness 3.2 mm) for 2.5 ml columns	20
<b>S1290</b>	Upper filter 90 µm for 2.5 ml columns	20
<b>S1310</b>	Upper filter 10 µm for 5 ml columns	20
<b>S1335</b>	Upper filter 35 µm (thickness 1.5 mm) for 5 ml columns	20
<b>S133532</b>	Upper filter 35 µm (thickness 3.2 mm) for 5 ml columns	20
<b>S1390</b>	Upper filter 90 µm for 5 ml columns	20
<b>S1410</b>	Upper filter 10 µm for 10 ml columns	20
<b>S143532</b>	Upper filter 35 µm (thickness 3.2 mm) for 10 ml columns	20
<b>S1490</b>	Upper filter 90 µm for 10 ml columns	20

## Column sets

Order No.	Description	Amount
<b>S10011</b>	Set: Mobicol "F" M1050; 2.5 ml; 5 ml and 10 ml columns (10 µm pore size filter)	3 each
<b>S1001</b>	Set: Mobicol "F" M1050; 2.5 ml; 5 ml and 10 ml columns (35 µm pore size filter)	3 each
<b>S10019</b>	Set: Mobicol "F" M1050; 2.5 ml; 5 ml and 10 ml columns (90 µm pore size filter)	3 each
<b>S10031</b>	Set upper filters 10 µm pore size for S10011	1 set
<b>S1003</b>	Set upper filters 35 µm pore size for S1001	1 set
<b>S10039</b>	Set upper filters 90 µm pore size for S10019	1 set






## Related Products: Matrices

A wide range of matrices, e.g. for the purification of His- or GST-tagged proteins, antibodies, nucleotide-binding proteins and phosphor-aminoacid-binding proteins, are available to be used in our Mobicol® or plastic columns.

An overview is listed in the table below. Please contact us for more information or visit our website.

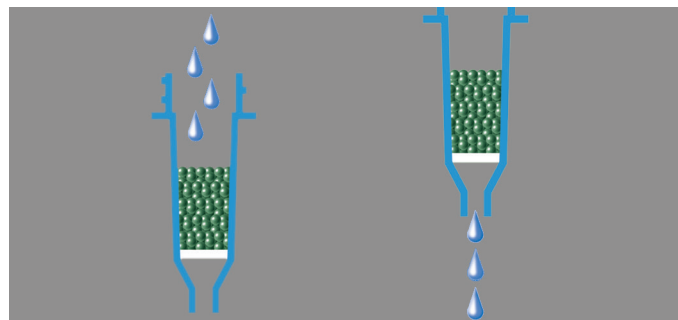

Order No.	Description	Amount
AC-201-10-JB	6% CL-Glutathione ChroMatrix™	10 ml
AC-303-25-JB	High Density Nickel Agarose	25 ml
AC-304-25-JB	Low Density Nickel Agarose	25 ml
AC-305-25-JB	High Density Zinc Agarose	25 ml
AC-202-10-JB	6% CL-Nickel ChroMatrix™	10 ml
AC-309-5-JB	Protein A Agarose	5 ml
AC-101S-JB	Immobilized gamma-Aminophenyl-ATP (C10-spacer), Agarose	1 ml
AC-102S-JB	Immobilized gamma-Aminophenyl-ATP (no spacer), Agarose	1 ml
AC-106S-JB	Immobilized gamma-Amino-octyl-GTP, Agarose	1 ml
AC-107S-JB	Immobilized gamma-Amino-octyl-UTP, Agarose	1 ml
AC-108S-JB	Immobilized gamma-Amino-octyl-CTP, Agarose	1 ml
AC-109S-JB	Immobilized gamma-Amino-octyl-ITP, Agarose	1 ml
AC-110S-JB	Immobilized gamma-Amino-octyl-XTP, Agarose	1 ml
AC-111S-JB	Immobilized gamma-Amino-octyl-dATP, Agarose	1 ml
AC-115S-JB	Immobilized gamma-Amino-octyl-dTTP, Agarose	1 ml
AC-119S-JB	Immobilized gamma-Amino-hexyl-CTP, Agarose	1 ml
AC-120S-JB	Immobilized gamma-Amino-hexyl-ITP, Agarose	1 ml
AC-121S-JB	Immobilized gamma-Amino-hexyl-XTP, Agarose	1 ml
AC-122S-JB	Immobilized gamma-Amino-hexyl-dATP, Agarose	1 ml
AC-124S-JB	Immobilized gamma-Amino-hexyl-dUTP, Agarose	1 ml
AC-137S-JB	Immobilized O-Phospho-L-Serine (C10-spacer), Agarose	1 ml
AC-139S-JB	Immobilized O-Phospho-L-Threonine (C10-spacer), Agarose	1 ml
AC-103S-JB	Immobilized O-Phospho-L-Tyrosine (C10-spacer), Agarose	1 ml
AC-104S-JB	Immobilized O-Phospho-L-Tyrosine (no spacer), Agarose	1 ml
AC-142S-JB	Immobilized 2'/3'-EDA-m7GTP, Agarose	1 ml
AC-144S-JB	Immobilized 2'/3'-EDA-ADP, Agarose	1 ml
AC-145S-JB	Immobilized N6-(6-Amino)hexyl-AMP, Agarose	1 ml
AC-146S-JB	Immobilized 8-Amino-hexyl-cAMP, Agarose	1 ml
AC-147S-JB	Immobilized 2'-EDA-cAMP, Agarose	1 ml
AC-148S-JB	Immobilized 2'-EDA-cGMP, Agarose	1 ml





## References

- Jinhi Ahn, Jason T. Wong, and Robert S. Molday The Effect of Lipid Environment and Retinoids on the ATPase Activity of ABCR, the Photoreceptor ABC Transporter Responsible for Stargardt Macular Dystrophy *J. Biol. Chem.*, Jun 2000; 275: 20399.
- Véronique Amiard, Annette Morvan-Bertrand, Jean-Pierre Billard, Claude Huault, Felix Keller, and Marie-Pascale Prud'homme Fructans, But Not the Sucrosyl-Galactosides, Raffinose and Loliose, Are Affected by Drought Stress in Perennial Ryegrass Plant *Physiology*, Aug 2003; 132: 2218.
- Pierre-Olivier Angrand, Inmaculada Segura, Pamela Völkel, Sonja Ghidelli, Rebecca Terry, Miro Brajenovic, Kristina Vintersten, Rüdiger Klein, Giulio Superti-Furga, Gerard Drewes, Bernhard Kuster, Tewis Bouwmeester, and Amparo Acker-Palmer Transgenic Mouse Proteomics Identifies New 14-3-3-associated Proteins Involved in Cytoskeletal Rearrangements and Cell Signaling *Mol. Cell. Proteomics*, Dec 2006; 5: 2211 - 2227.
- Denis Bayle, Sabine Wängler, Thomas Weitzenegger, Wolfram Steinhilber, Jürgen Volz, Michael Przybylski, Klaus P. Schäfer, George Sachs, and Klaus Melchers Properties of the P-Type ATPases Encoded by the copAP Operons of *Helicobacter pylori* and *Helicobacter felis* *J. Bacteriol.*, Jan 1998; 180: 317 - 329.
- Miro Brajenovic, Gerard Joberty, Bernhard Kuster, Tewis Bouwmeester, and Gerard Drewes Comprehensive Proteomic Analysis of Human Par Protein Complexes Reveals an Interconnected Protein Network *J. Biol. Chem.*, Mar 2004; 279: 12804 - 12811.
- Jan Brix, Klaus Dietmeier, and Nikolaus Pfanner Differential Recognition of Preproteins by the Purified Cytosolic Domains of the Mitochondrial Import Receptors Tom20, Tom22, and Tom70 *J. Biol. Chem.*, Aug 1997; 272: 20730.
- Jan Brix, Stefan Rüdiger, Bernd Bukau, Jens Schneider-Mergener, and Nikolaus Pfanner Distribution of Binding Sequences for the Mitochondrial Import Receptors Tom20, Tom22, and Tom70 in a Presequence-carrying Preprotein and a Non-cleavable Preprotein *J. Biol. Chem.*, Jun 1999; 274: 16522.
- Alejandro Ferrando, Zsuzsana Koncz-Kálmán, Rosa Farràs, Antonio Tiburcio, Jeff Schell, and Csaba Koncz Detection of in vivo protein interactions between Snf1-related kinase subunits with intron-tagged epitope-labelling in plants cells *Nucleic Acids Res.*, Sep 2001; 29: 3685 - 3693.
- J. Fleck Chemical Fate of a Metamorphic Inducer in Larvae-like Buds of the Cnidarian *Cassiopea andromeda* *Biol. Bull.*, Feb 1998; 194: 83.
- Albrecht Gruhler, Jesper V. Olsen, Shabaz Mohammed, Peter Mortensen, Nils J. Færgeman, Matthias Mann, and Ole N. Jensen Quantitative Phosphoproteomics Applied to the Yeast Pheromone Signaling Pathway *Mol. Cell. Proteomics*, Mar 2005; 4: 310 - 327.
- Klaus Hellmuth, Denise M. Lau, F. Ralf Bischoff, Markus Künzler, Ed Hurt, and George Simos. Cell and Organelle Structure and Assembly: Yeast Los1p has Properties of an Exportin-Like Nucleocytoplasmic Transport Factor for tRNA.
- Munir Iqbal, Emma Poole, Stephen Goodbourn, and John W. McCauley Role for Bovine Viral Diarrhea Virus Erns Glycoprotein in the Control of Activation of Beta Interferon by Double-Stranded RNA *J. Virol.*, Jan 2004; 78: 136 - 145.
- Irvine SQ Whole-mount in situ hybridization of small invertebrate embryos using laboratory mini-columns *Biotechniques*, Dec 2007; 43(6): 764, 766, 768.
- Florian Krötz, Barbara Engelbrecht, Martin A. Buerkle, Florian Bassermann, Hanna Bridell, Torsten Gloe, Justus Duyster, Ulrich Pohl, and Hae-Young Sohn The Tyrosine Phosphatase, SHP-1, Is a Negative Regulator of Endothelial Superoxide Formation *J. Am. Coll. Cardiol.*, May 2005; 45: 1700 - 1706.
- Jelle Van Leene, Hilde Stals, Dominique Eeckhout, Geert Persiau, Eveline Van De Slijke, Gert, Van Isterdael, Annelies De Clercq, Eric Bonnet, Kris Laukens, Noor Remmerie, Kim Hendrickx, Thomas De Vijlder, Azmi Abdelkrim, Anne Pharazyn, Harry Van Onckelen, Dirk Inzé, Erwin Witters, and Geert De Jaeger A Tandem Affinity Purification-based Technology Platform to Study the Cell Cycle Interactome in *Arabidopsis thaliana* *Mol. Cell. Proteomics*, Jul 2007; 6: 1226 - 1238.
- Jérémy Lothier, Bertrand Lasseur, Katrien Le Roy, André Van Laere, Marie-Pascale Prud'homme, Philippe Barre, Wim Van den Ende, and Annette Morvan-Bertrand Cloning, gene mapping, and functional analysis of a fructan 1-exohydrolase (1-FEH) from *Lolium perenne* implicated in fructan synthesis rather than in fructan mobilization *J. Exp. Bot.*, Jun 2007; 58: 1969 - 1983.
- A. Magrelli, K Langenkemper, C Dehio, J Schell, and A Spena Splicing of the *rolA* transcript of *Agrobacterium rhizogenes* in *Arabidopsis* *Science*, Dec 1994; 266: 1986 - 1988.
- Heidmarie Müllner, Günter Deutsch, Erich Leitner, Elisabeth Ingolic, and Günther Daum YEH2/YLR020c Encodes a Novel Steryl Ester Hydrolase of the Yeast *Saccharomyces cerevisiae* *J. Biol. Chem.*, Apr 2005; 280: 13321 - 13328.
- Karsten Niederhoff, Nadja M. Meindl-Beinker, Daniela Kerksen, Uta Perband, Antje Schäfer, Wolfgang Schliebs, and Wolf-H. Kunau Yeast Pex14p Possesses Two Functionally Distinct Pex5p and One Pex7p Binding Sites *J. Biol. Chem.*, Oct 2005; 280: 35571 - 35578.
- Stephenie Paine-Saunders, Beth L. Viviano, Aris N. Economides, and Scott Saunders Heparan Sulfate Proteoglycans Retain Noggin at the Cell Surface. A Potential Mechanism for Shaping Bone Morphogenetic Protein Gradients *J. Biol. Chem.*, Jan 2002; 277: 2089 - 2096.
- K. Pfeifer, M Bachmann, HC Schroder, BE Weiler, D Ugarkovic, T Okamoto, and WE Muller Formation of a small ribonucleoprotein particle between Tat protein and trans-acting response element in human immunodeficiency virus-infected cells *J. Biol. Chem.*, Aug 1991; 266: 14620 - 14626.
- Werner G. Purschke, Falko Radtke, Frank Kleinjung, and Sven Klussmann A DNA Spiegelmer to staphylococcal enterotoxin B *Nucleic Acids Res.*, Jun 2003; 31: 3027 - 3032.
- Juri Rappsilber, Paul Ajuh, Angus I. Lamond, and Matthias Mann SPF30 Is an Essential Human Splicing Factor Required for Assembly of the U4/U5/U6 Tri-small Nuclear Ribonucleoprotein into the Spliceosome *J. Biol. Chem.*, Aug 2001; 276: 31142 - 31150.
- B. Ashok Reddy, Prashanth Kumar Bajpe, Andrew Bassett, Yuri M. Moshkin, Elena Kozhevnikova, Karel Bezstarosti, Jeroen A.A. Demmers, Andrew A. Travers, and C. Peter Verrijzer Drosophila Transcription Factor Tramtrack69 Binds MEP1 to Recruit the Chromatin Remodeler NuRD *Mol. Cell. Biol.*, Aug 2010; 30: 1128/ MCB.00266-10.



## Chapter II

**Spin Columns for your  
own Matrices  
or other Compounds**



## MobiSpin Column “F”

### Advantages

**Spin column chromatography offers many advantages over traditional liquid chromatography:**



- Comes empty - fill in your own material (matrix, gel slice etc.)
- Easy handling: load sample, spin and collect the purified product
- No sample dilution
- Reproducible results with fast protocols
- One application in less than 4 minutes
- Numerous samples can be processed simultaneously
- Large number of applications
- Compatible with laboratory standard

### Background


The MobiSpin column “F” is designed for a wide variety of applications for nucleic acid and protein extraction and purification. When choosing the MobiSpin column “F” for a particular application, a suitable resin must be selected and filled into the mini-column. For an appropriate resin selection, impurity

versus sample size and the anticipated yield must be considered.

The MobiSpin column “F” comes with inserted small or large filter of 10 µm pore size, screw cap and snap-off plug. It can be used without matrix (e.g. to extract trypsinated proteins out of a gel slice for Mass Spec analysis) or filled with a matrix of your choice (see also Chapter III).

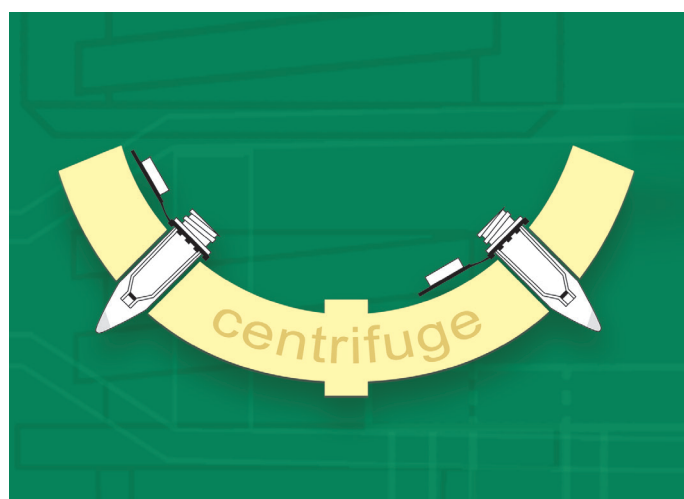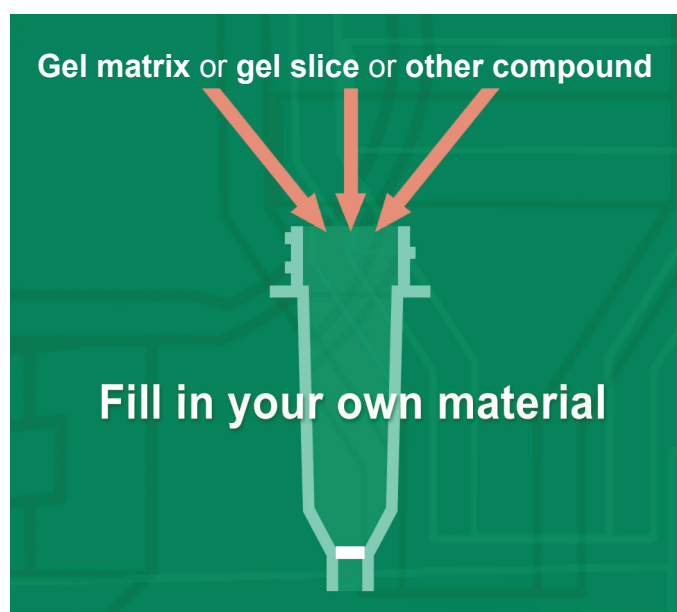

MobiSpin Column “F” can be centrifuged in 1.5 ml or 2 ml tubes

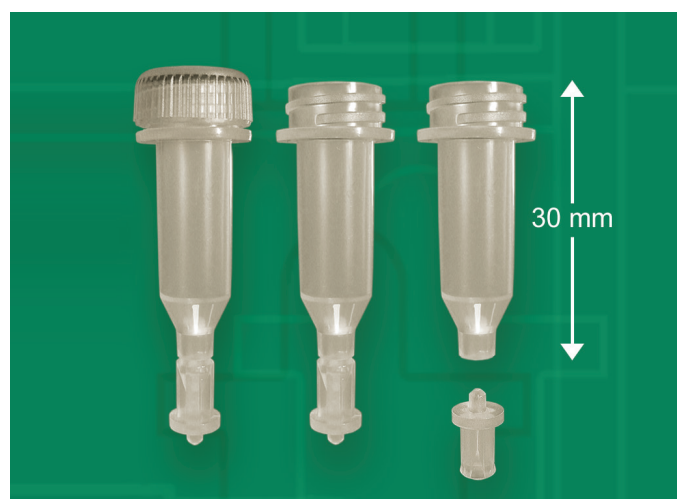


The MobiSpin column is autoclavable at 110 °C for 10 min. Close the column with the tightly closing screw cap. Snap off the outlet plug and use its reverse side for closing.



## Special Application of MobiSpin Column “F”

### Preparation of trypsinated protein fragments for Mass Spec

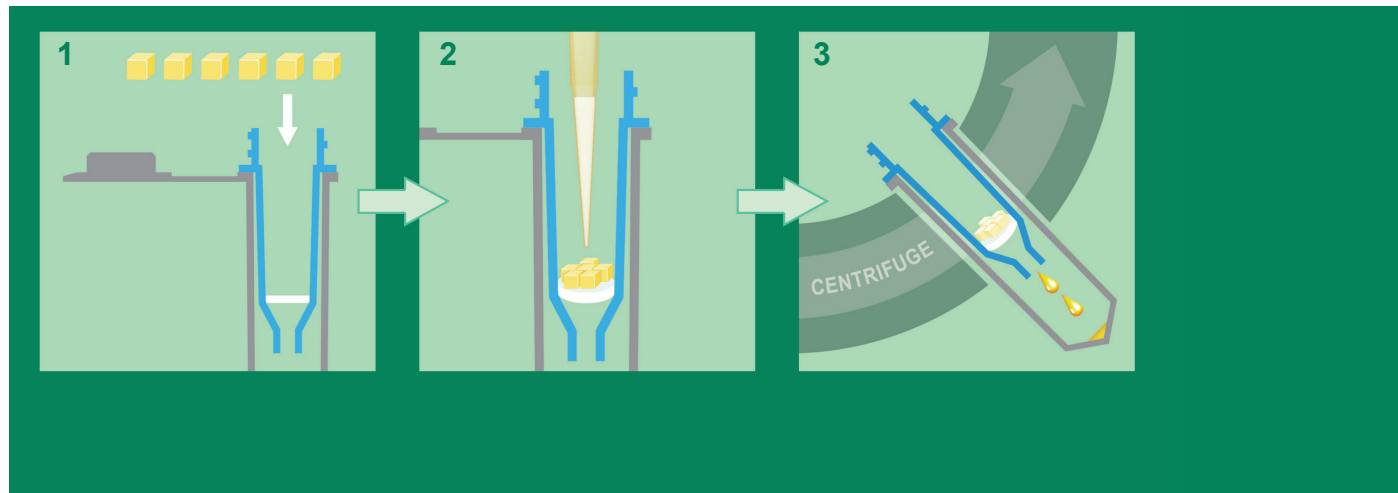

For mass spectrometry analysis, proteins are separated by polyacrylamide gel electrophoresis (PAGE) and a single gel slice containing a protein of interest is cut out of the gel. The slice is cut into pieces of approximately 2 x 2 mm and then placed onto the inserted filter of a MobiSpin column “F” (M105210S) (1). After several washing and drying steps involving centrifugation, trypsin is added to the gel pieces inside of the column, diffuses into them and digests the protein overnight at 37 °C into small fragments (2). The protein fragments are eluted by a few spinning

steps (3). After drying the collected spin elutes can be analyzed by mass spectrometry. The complete protocol is available on our web page ([www.mobitec.com](http://www.mobitec.com)). The peptide mass fingerprints provided by this protocol are consistent with results obtained by other common mass spec sample preparations.

This protocol will be suitable for proteases other than trypsin as well.

# Ask for Free Samples!

### MobiSpin Column “F”

Order No.	Description	Amount
<b>M105010S</b>	MobiSpin Column “F” with fixed outlet plug, inserted small 10 µm filter and screw cap	50
<b>M105210S</b>	MobiSpin Column “F” with fixed outlet plug, inserted large 10 µm filter and screw cap	50





## References

Demian Cazalla, Kathryn Newton, and Javier F. Cáceres A Novel SR-Related Protein Is Required for the Second Step of Pre-mRNA Splicing *Mol. Cell. Biol.*, Apr 2005; 25: 2969 - 2980.

B. K. Choi, B J Paster, F E Dewhirst, and U. B. Göbel Diversity of cultivable and uncultivable oral spirochetes from a patient with severe destructive periodontitis. *Infect. Immun.*, May 1994; 62: 1889 - 1895.

Holger Eubel, Etienne H. Meyer, Nicolas L. Taylor, John D. Busell, Nicholas O'Toole, Joshua L. Heazlewood, Ian Castleden, Ian D. Small, Steven M. Smith, and A. Harvey Millar Novel Proteins, Putative Membrane Transporters, and an Integrated Metabolic Network Are Revealed by Quantitative Proteomic Analysis of Arabidopsis Cell Culture Peroxisomes *Plant Physiology*, Dec 2008; 148: 1809 - 1829.

J. Hedegaard, SA Steffensen, N Nørskov-Lauritsen, KK Mortensen, and HU Sperling-Petersen Identification of Enterobacteriaceae by partial sequencing of the gene encoding translation initiation factor 2 *Int. J. Syst. Bacteriol.*, Oct 1999; 49: 1531.

Yoko Kita, Yoshiaki Miura, Jun-ichi Furukawa, Mika Nakano, Yasuhiro Shinohara, Masahiro Ohno, Akio Takimoto, and Shin-Ichiro Nishimura Quantitative Glycomics of Human Whole Serum Glycoproteins Based on the Standardized Protocol for Liberating N-Glycans *Mol. Cell. Proteomics*, Aug 2007; 6: 1437 - 1445.

Angeles Ortega, Martina Niksic, Angela Bachi, Matthias Wilm, Lucas Sánchez, Nicholas Hastie, and Juan Valcárcel Biochemical Function of Female-Lethal (2)D/Wilms' Tumor Suppressor-1-associated Proteins in Alternative Pre-mRNA Splicing, *J. Biol. Chem.*, Jan 2003; 278: 3040 - 3047.

Shaun Peters, Sagadevan G. Mundree, Jennifer A. Thomson, Jill M. Farrant, and Felix Keller Protection mechanisms in the resurrection plant *Xerophyta viscosa* (Baker): both sucrose and raffinose family oligosaccharides (RFOs) accumulate in leaves in response to water deficit *J. Exp. Bot.*, Jun 2007; 58: 1947 - 1956.

Stephanie Pfander, Roberto Fiammengo, Srečko I. Kirin, Nils Metzler-Nolte, and Andres Jäschke Reversible site-specific tagging of enzymatically synthesized RNAs using aldehyde-hydrazine chemistry and protease-cleavable linkers *Nucleic Acids Res.*, Feb 2007; 35: e25.

Sonja Reiland, Gaëlle Messerli, Katja Baerenfaller, Bertran Gerlits, Anne Endler, Jonas Grossmann, Wilhelm Grissem, and Sacha Baginsky Large-Scale Arabidopsis Phosphoproteome Profiling Reveals Novel Chloroplast Kinase Substrates and Phosphorylation Networks *Plant Physiology*, Jun 2009; 150: 889 - 903.

Symeon Siniouoglou, Ed C. Hurt, and Hugh R. B. Pelham Psr1p/Psr2p, Two Plasma Membrane Phosphatases with an Essential DXDX(T/V) Motif Required for Sodium Stress Response in Yeast *J. Biol. Chem.*, Jun 2000; 275: 19352.

W. Voos, B. D. Gambill, S. Laloraya, D. Ang, E. A. Craig, and N. Pfanner Mitochondrial GrpE is present in a complex with hsp70 and preproteins in transit across membranes. *Mol. Cell. Biol.*, Oct 1994; 14: 6627 - 6634.

F. T. Vreede and G. G. Brownlee Influenza Virion-Derived Viral Ribonucleoproteins Synthesize both mRNA and cRNA *In Vitro J. Virol.*, Mar 2007; 81: 2196 - 2204.

Juan Zou, Lorna Henderson, Vicky Thomas, Patricia Swan, A. Neil Turner, and Richard G. Phelps Presentation of the Good-pasture Autoantigen Requires Proteolytic Unlocking Steps That Destroy Prominent T Cell Epitopes *J. Am. Soc. Nephrol.*, Mar 2007; 18: 771 - 779.





# Chapter III


## **Matrix-filled MobiSpin Columns for Purification of Nucleic Acids**

## Pre-packed MobiSpin Columns for Separation of DNA and RNA from Small Contaminants

Removing salts and other small molecules from biomolecules is often essential for downstream processes. MobiSpin Columns are developed for rapid and efficient routine separation tasks regarding DNA and RNA.

### Features

- Compatibility with **laboratory standard**
- Columns are **pre-packed** and **ready to use**
- **Very fast** procedure: Sample loading and recovery in less than 3 minutes
- **Easy handling:**
  1. Equilibrate with MobiSpin or Storage Buffer, spin, and discard the flow through
  2. Load sample, spin, and collect the purified product
- Numerous samples can be **processed simultaneously**
- **Excellent results:** Exceptional high DNA/RNA recovery with > 95% desalting efficiency
- **No sample dilution**

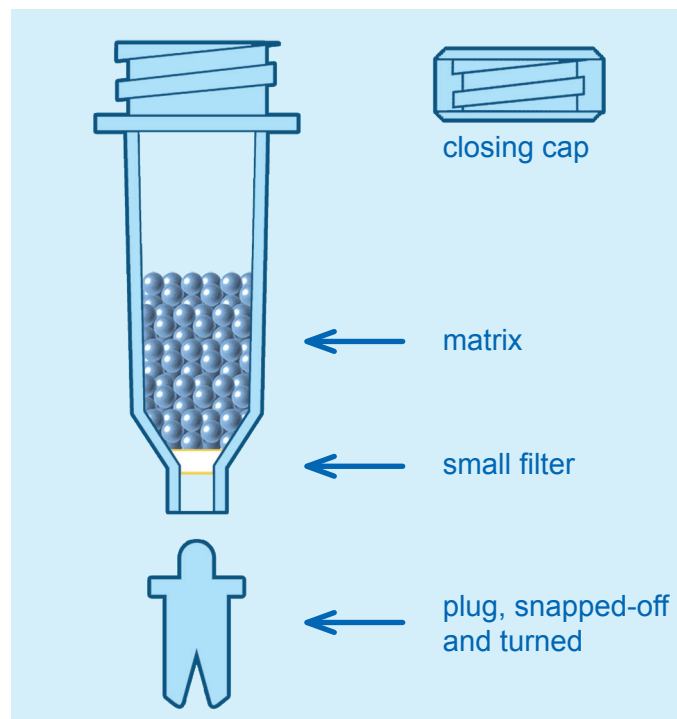

### Pre-packed MobiSpin Column

The pre-packed MobiSpin Column is closed with a screw cap and contains a matrix comprised of microscopic beads of a porous material.

The tip of the MobiSpin Column is a fixed snap-off plug. Prior to the first centrifugation step it has to be bent down and removed. The plug can be reused for closing by turning it upside down.








### Functional Principle: Size-Exclusion Chromatography

MobiSpin Columns, pre-packed with diverse matrices, are designed for a wide variety of separation tasks including nucleic acid purification. They function by the proven principle of size-exclusion chromatography. The columns combine the effectiveness of gel filtration with speed of centrifugation.

The pore size of the filled-in matrix determines which molecules are small enough to enter the pores of the matrix beads and which molecules are too large. In this way nucleic acid molecules larger than the pore size are excluded from the resin and located in the

void volume. These molecules move quickly through the matrix bed when the column is centrifuged. Molecules smaller than the pore size, like hydrated salt ions, enter the pores of the matrix beads and are held back.

A detailed overview on the functional principle is given in the figure below.

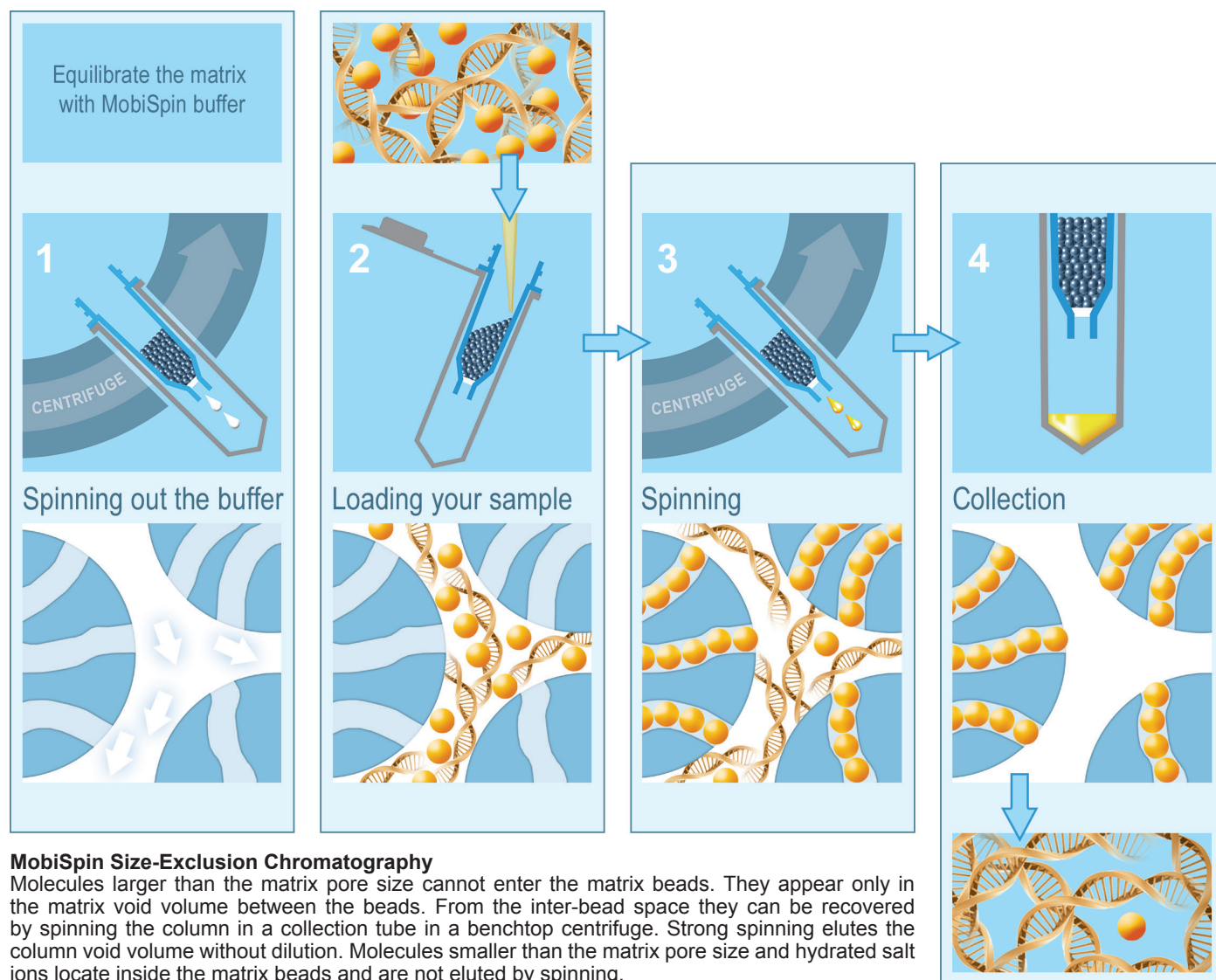

## Applications

Pre-packed MobiSpin Columns for nucleic acid purification are available with three different matrices.


- Silica: Desalting MobiSpin Columns
- Sephacryl®: MobiSpin S-Columns
- Sephadex®: MobiSpin G-Columns

The resins have differing suitability regarding their usage. We recommend these pre-packed MobiSpin columns as rapid and efficient tool for the following applications.

### Desalting and Buffer Exchange

**Desalting MobiSpin Columns** quickly remove small molecules like salt ions from DNA/RNA (even from dsDNA fragments <25 bp and from dNTPs) with good recovery efficiency. The columns are pre-packed with a novel silica resin comprised of uniform microscopic beads of a porous, form-stable material. Since this new matrix material works without moisture expansion, resins are allowed to fall dry without any loss in purification efficiency. This is in contrast to Sephacryl®, Sephadex®, and Sepharose® resins which swell in buffer. Therefore, Desalting MobiSpin Columns give reproducible results of high quality and are much easier to handle. Furthermore, the columns are long-term storable when kept dry.

As an alternative, **MobiSpin G-Columns** pre-packed with conventional Sephadex® G-50 matrix may be used. This hydrophilic matrix swells in buffer and acquires its chromatographic properties only under wet conditions. In this state beads have a pore size of 700 Da that allow hydrated salt ions to enter into the pores, while DNA/RNA > 20 bases and most other biomolecules stay outside.

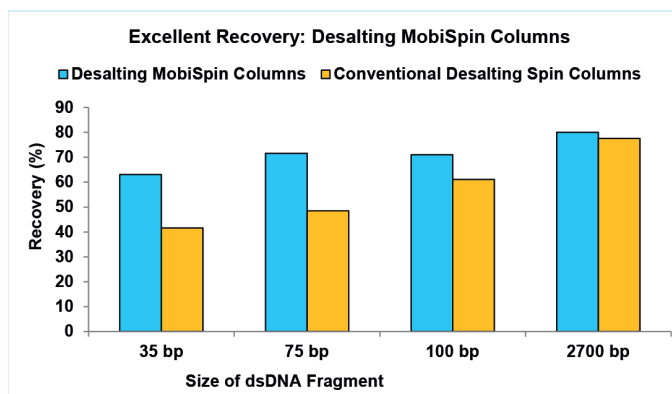

### DNA Recovery of Desalting MobiSpin Columns and Conventional Desalting Spin Columns

In all cases, a high desalting efficiency of > 90% could be achieved. The recovery of large-sized fragments > 2kb was nearly comparable (about 80%) with both kinds of spin columns. But with small-sized fragments ≤ 100 bp the Desalting MobiSpin Columns showed a considerably improved recovery.

### Fast Removal of Dye Terminators or Unincorporated Labeled Nucleotides

**MobiSpin S-Columns** pre-packed with Sephacryl® matrices (S-200, S-300, S-400) or **MobiSpin G-Columns** pre-packed with Sephadex® G-50 matrix are particularly suited for the removal of fluorescent dye dideoxyterminators (e.g., Cy5/Cy3 nucleotides) from cycle sequencing reactions.

Furthermore these columns are convenient for the extraction of unincorporated labeled nucleotides (dye-labeled or radiolabeled dNTPs or ddNTPs) from DNA labeling reactions, e.g., PCR probe labeling, Nick Translation, or DNA end-labeling. The purified DNA is applicable to downstream applications like FISH (fluorescence *in situ* hybridization) or Southern/Northern blotting. The removal of unincorporated labeled nucleotides is a precondition for determining the DNA labeling rate. For good recovery rates, labeled DNA fragments must be at least 20 bp in length.

**Pre-packed and  
ready-to-use**

Sephacryl®, Sephadex® and Sepharose® are registered trademarks of GE Healthcare.



## Applications

### Removal of dNTPs, Oligos, and Salt

**MobiSpin S-Columns** pre-packed with Sephacryl® matrices (S-200, S-300, S-400) are particularly suited for purification of PCR and other enzymatic DNA reactions. With these matrices nucleotides, oligonucleotides, and buffer (but not enzymes) will be removed. The larger the pore size of the resin (S 400 > S-300 > S-200), the greater the purity and lower the yield of the product.

For obtaining optimal results please consider: the smallest product being purified should be at least 20 times larger than the largest impurity.

### Note:

For choosing the appropriate Sephacryl® matrix, the exact size of the nucleic acid fragment to be purified and the sample volume has to be considered. As a guideline we recommend using MoBiTecs "Sample Volume Guide" that is given in the **MobiSpin S-Columns Handbook** which is available on our website.

### Pre-packed MobiSpin Columns - Overview

Column	Matrices	Applications
Desalting MobiSpin	Silica	Desalting and buffer exchange especially of small sized nucleic fragments and dNTPs
MobiSpin S (S-200; S-300; S-400)	Sephacryl®	Buffer exchange, nucleotide removal, plasmid, oligo and PCR reaction purification, removal of unincorporated dyes or dye terminators
MobiSpin G (G-50)	Sephadex®	Removal of unincorporated dyes or dye terminators, desalting and buffer exchange

**Sample loading and  
recovery in less  
than 3 minutes**



## Order Information, Shipping and Storage

### Silica-Based Columns

Order No.	Description	Amount
SCO100	Desalting MobiSpin	20 Columns
SCO110	Desalting MobiSpin	100 Columns

Including MobiSpin Buffer (10 mM Tris/HCl pH 7.6; 1 mM EDTA) for equilibration; shipped at RT; store columns at RT and buffer at 4 °C

### Sephacryl®-Based Columns

Order No.	Description	Amount
SCO200	MobiSpin S-200	20 Columns
SCO210	MobiSpin S-200	100 Columns
SCO234	MobiSpin S-200, S-300, S-400	3x 10 Columns
SCO300	MobiSpin S-300	20 Columns
SCO310	MobiSpin S-300	100 Columns
SCO400	MobiSpin S-400	20 Columns
SCO410	MobiSpin S-400	100 Columns

Equilibrated in MobiSpin Buffer (10 mM Tris/HCl pH 7.6; 1 mM EDTA); shipped at RT; store columns at 4 °C

### Sephadex®-Based Columns

Order No.	Description	Amount
SCO500	MobiSpin G-50	20 Columns
SCO510	MobiSpin G-50	100 Columns

Equilibrated in MobiSpin Buffer (10 mM Tris/HCl pH 7.6; 1 mM EDTA) shipped at RT; store columns at 4 °C

## References

Frei A. P. *et al.* (2013) Ligand-based receptor identification on living cells and tissues using TRICEPS; *Nature Protocols* 8, 1321–1336

Maeda K. *et al.* (2014) A generic protocol for the purification and characterization of water-soluble complexes of affinity-tagged proteins and lipids; *Nature Protocols* 9, 2256–2266

Marini F. and Wood R. D. (2002) A human DNA helicase homologous to the DNA cross-link sensitivity protein Mus308; *J. Biol. Chem.* 277, 8716–8723

Porubsky St. *et al.* (2004) Influence of native and hypochlorite-modified low density lipoprotein on gene expression in human proximal tubular epithelium; *Am. J. Pathol.* 164, 2175–2187

Sachsenröder J. *et al.* (2014) The general composition of the faecalvirome of pigs depends on age, but not on feeding with a probiotic bacterium; *PLOS ONE* 9 (2) e88888

Schneider J. *et al.* (2004) Systematic analysis of T7 RNA polymerase based *in vitro* linear RNA amplification for use in microarray experiments; *BMC Genomics* 5 (29) 1471–2164

Trojnar E. *et al.* (2009) The first complete genome sequence of a chicken group A rotavirus indicates independent evolution of mammalian and avian strains; *Virology* 386 (2) 325–333



# Chapter IV

## **Mobicols for Affinity Purification of Biomolecules**

### **Fractionation and Purification of Antibodies**

Purification of antibodies using immobilized Protein A / Protein G

### **Immobilized Metal Ion Affinity Chromatography (IMAC)**

Purification of recombinant polyhistidine-tagged proteins  
using pre-packed Ni-IDA columns





## Protein A and G columns

### Advantages

- Time saving: antibody purification in less than 20 minutes
- High-quality controlled glass porous matrix is more rigid and durable than conventional matrices
- Starter kits with required buffers, columns and accessories available
- The Protein A and G columns are ideally suited for the efficient and rapid purification of antibodies with minimal pre-treatment
- Protein A and G columns have good flow properties
- Capacity of these Protein A and G columns for murine IgG at physiological pH and salt concentration are greater than that for other Protein A and G adsorbents

### Background

The protein is immobilized on controlled pore glass which is more rigid and durable than conventional polymeric (for example agarose) matrices. This glass also has a very uniform internal pore size distribution which falls within a very narrow range, resulting in uniform rates of diffusion and hence rapid transfer between solution and solid phase. The protein molecules are orientated and immobilized on the glass surface such that pairs of Protein A and G molecules can interact with the binding regions of each of the two Fc portions of the immunoglobulin at the same time. This results in very high binding capacities, especially for those immunoglobulins which typically do not bind well to conventional immobilized Protein A and G and where this bivalent binding is an essential prerequisite (for example mouse IgG<sub>1</sub>). The increased capacity does not, however, result in a higher association constant ( $K_a$ ) and therefore, elution conditions remain unaltered.



**Columns with immobilized Protein A and G may be used for many applications including:**

- Purification of monoclonal antibodies
- Selective purification of antibody fragments
- Purification or removal of polyclonal IgG from serum
- Separation of IgG sub-classes using a stepwise pH gradient

The plastic columns with immobilized Protein A and G facilitate highly efficient, cost effective purification for both monoclonal and polyclonal antibodies. These columns bind antibodies under binding buffer conditions, the antibodies are eluted from the column under elution buffer conditions, and the columns are then washed and stored under storage buffer. They can be stored and reused many times.


**The immobilization of Protein A and G has been developed to satisfy four critical factors:**

- Orientation of the ligand
- Distribution of the ligand
- Stability of the immobilized ligand
- Elimination of non-specific surface interactions

Immobilized onto a solid glass support and packed into columns, Protein A and G purifies total IgG from crude protein mixtures such as serum or ascites fluid.

## Antibody purification by Protein A and G columns


Polyclonal antibodies in plasma or serum, or monoclonal antibodies in cell culture supernatant are loaded onto the Protein A or G columns. There, the Fc domains of the antibodies are captured by the matrix immobilized Protein A or G (Loading) while other proteins are not bound. All contaminants other than the antibodies are washed off (Washing) removing unbound or weakly bound material from the column. Then, the buffer conditions are changed to lower pH values (down to about pH 3) so that Protein A or G cannot bind the Fc antibody domains anymore and the purified antibodies are eluted (Elution).

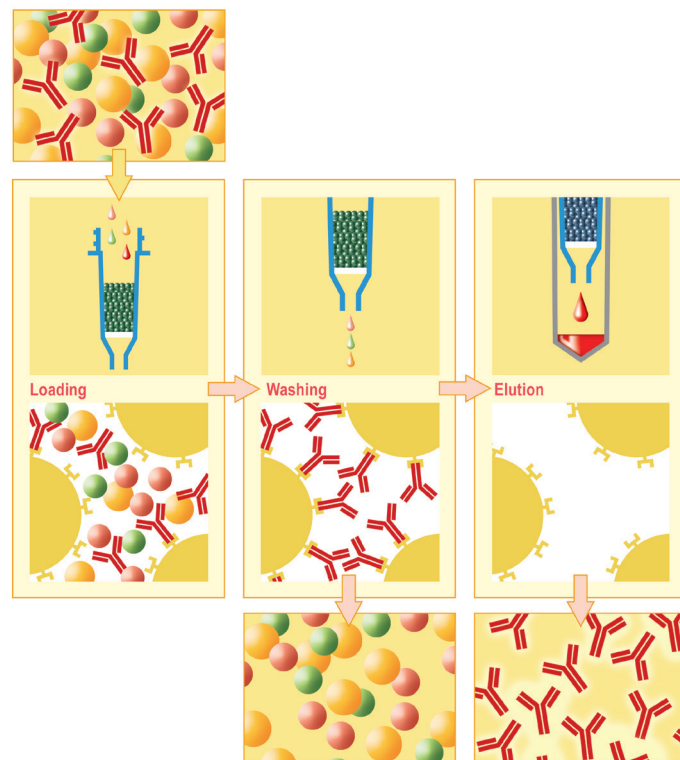

## Example

### Purification of polyclonal antibodies from rabbit serum

The elution yield is determined by SDS-PAGE gel analysis of the washing and elution steps

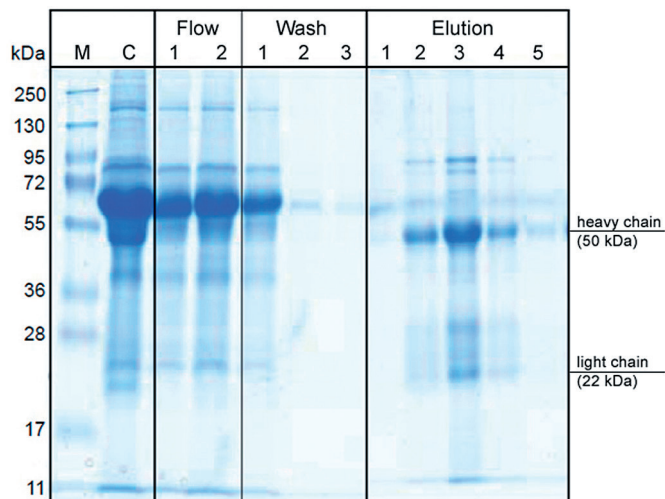

SDS-PAGE analysis of washing and elution steps. M: protein markers, C: loading, Flow: flow through 1 and 2, Wash: washing steps 1, 2 and 3, Elution: elution steps 1, 2, 3, 4 and 5. 1:4 diluted rabbit serum was loaded onto the column (C) and the flow through analyzed (Flow). For washing (Wash), about 5 ml PBS buffer (10 column volumes) was loaded onto the column and the flow through was collected in 1.5 ml aliquots; 3 aliquots were loaded on the gel. For elution (Elution), 5 ml elution buffer was loaded onto the column and the flow through was collected in 0.5 ml aliquots; the first 5 aliquots were loaded on the gel. The large amount of the dominant 65 kDa protein runs through the column and does not bind. The antibody heavy chain (50 kDa) does bind to the column (band present in lane C but not in Flow and Wash lanes). The antibody heavy (50 kDa) and light (22 kDa) chains elute in elution aliquots 2 to 4 and mainly in elution step 3.





## Applications

### Protein A, Protein G

Protein A and Protein G are bacterial proteins with high specific binding to the Fc (non-antigen binding) domain of many classes of antibodies. Protein A and G can be used to affinity purify immunoglobulins (primarily IgGs) or for the separation of Fc from Fab fragments. The binding characteristics of Protein A and G differ:

Protein A has a broad species reactivity. It binds well to IgG from human, rabbit, cow and guinea pig but weak to monoclonal antibodies.

Protein G has a broader species reactivity. It shows much stronger binding to IgG from mouse, rat and goat and also stronger binding to monoclonal antibodies.

Antibodies are eluted from the Protein A and G columns by low pH elution buffer.

### Features

- Purification of monoclonal antibodies
- Selective purification of antibody fragments
- Purification or removal of polyclonal IgG from serum
- Separation of IgG sub-classes using a stepwise pH gradient
- Protein A for IgG from human (IgM, IgA, IgE, IgG<sub>1</sub>, IgG<sub>2</sub>, IgG<sub>4</sub>), cat, dog, mouse, rabbit, pig, Guinea pig
- Protein G for IgG from human (IgG<sub>1</sub>, IgG<sub>2</sub>, IgG<sub>3</sub>, IgG<sub>4</sub>), cow, goat, horse, mouse (not IgM), pig, rabbit, rat and sheep

### Order Information, Shipping & Storage

Order No.	Description	Amount
<b>PR-PAK002</b>	Protein A Starter Kit, 2 columns, buffers, 2 ultrafiltration spin columns and tubes	Kit
<b>PR-PAK005</b>	Protein A Kit with 5 columns	Kit
<b>PR-PAK010</b>	Protein A Kit with 10 columns	Kit
<b>PR-PGK002</b>	Protein G Starter Kit, 2 columns, buffers, 2 ultrafiltration spin columns and tubes	Kit
<b>PR-PGK005</b>	Protein G Kit with 5 columns	Kit
<b>PR-PGK010</b>	Protein G Kit with 10 columns	Kit
<b>PR-PK001</b>	Buffers for Protein A/G Kits for 10 applications with ultrafiltration spin filters and tubes	Kit


Shipped at RT; store at 4 °C

Protein	A	G
Human IgG <sub>1</sub>	++	++
Human IgG <sub>2</sub>	++	++
Human IgG <sub>3</sub>	-	++
Human IgG <sub>4</sub>	++	++
Human IgA	+	-
Human IgD	+	-
Human IgE	+	-
Human IgM	+	-
Mouse IgG <sub>1</sub>	+	+
Mouse IgG <sub>2a</sub>	++	++
Mouse IgG <sub>2b</sub>	++	++
Mouse IgG <sub>3</sub>	+	++
Mouse IgM	+/-	-
Rat IgG	++	++
Rat IgG1	+/-	+
Rat IgG2a	+/-	++
Rat IgG2b	+/-	+
Rat IgG2c	+/-	+
Rat IgM	+/-	-
Rabbit IgG	++	++
Hamster IgG	+	++
Guinea pig IgG	++	+
Bovine IgG	+	+
Sheep IgG	+/-	+
Goat IgG	+/-	+
Pig IgG	++	++
Chicken IgG	-	+/-
<b>Fragments</b>		
Human Fab	+	+
Human F(ab') <sub>2</sub>	+	+
Human scFv	+	-
Human Fc	+	+
Human κ	-	-
Human λ	-	-
++ = strong affinity		
+ = Moderate/slight affinity		
+/- = Requires evaluation		
- = No affinity		

## Ready-to-Use Ni-IDA Columns – First Choice for Higher Yields of His- tagged Proteins

Ni-IDA Columns provide a fast and convenient routine tool for purification of recombinant polyhistidine-tagged proteins.

The form-stable silica matrix is precharged with  $\text{Ni}^{2+}$  ions and allows purification on the principle of Immobilized Metal Ion Affinity Chromatography (IMAC). Binding of proteins is based on the interaction between the polyhistidine tag of the recombinant protein and immobilized  $\text{Ni}^{2+}$  ions.

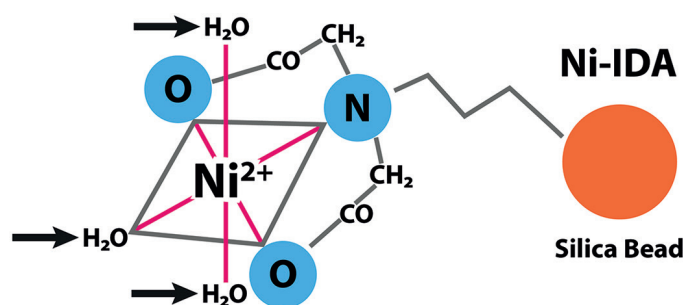


### Silica-Based Ni-IDA Matrix

The Ni-IDA silica is an affinity chromatography matrix for purifying recombinant proteins carrying a His-tag. Histidine residues in the His-tag bind to the vacant positions in the coordination sphere of the immobilized nickel ions which enables strong and efficient binding of target protein. In free Ni-IDA the vacant positions are usually occupied by water molecules (as shown above).

### Features

- Easy, fast, and cost-effective routine purification
- High binding capacity and high affinity
- Highest flexibility of purification:
  1. under native and denaturing conditions
  2. by gravity flow or by spinning
- Columns are **pre-packed** and **ready-to-use**
- Starting from diverse expression systems, e.g., *E. coli*, yeast, insect, and mammalian cells
- Universal use: suitable for small proteins, large protein complexes, proteins with low expression rates
- Improved target specificity by optimized silica-based  $\text{Ni}^{2+}$ -IDA matrix
- Imidazol-free loading and washing buffer
- Straightforward test of experimental protocols using His<sub>6</sub>-GFP bacterial expression vector

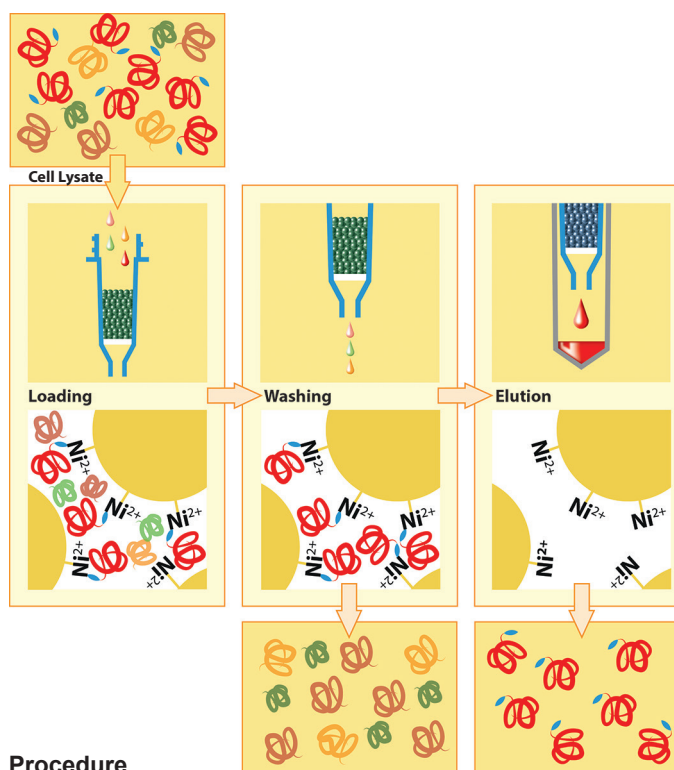

### Procedure

The purification of His-tagged proteins consists of 4 steps: cell lysis, loading & binding, washing, and elution. Cleared cell lysates are loaded onto the matrices. His-tagged proteins are bound, and other proteins pass through the matrix. After washing, His-tagged proteins are eluted with buffer from the matrix.

**Simply replace your current  
Ni-NTA products,  
no optimization or protocol  
change necessary!**

## MobiSpin Ni-IDA Columns

### Features

- For protein purification by spinning
- Very high binding capacity: up to 12 mg protein per spin column
- Excellent protein recovery rate of > 90%
- Simultaneous processing of multiple samples



## MoBiTec Ni-IDA Columns

### Features

- For protein purification by gravity flow
- Maximal binding capacity: up to 90 mg protein per column
- Protein recovery rate of > 80%
- Columns are long-term storable and reusable



### Background

The chelating group of the Ni-IDA resin is based on IDA (iminodiacetic acid), which enables strong and efficient binding of target protein onto the IMAC matrix.

In contrast to traditional IDA matrices, MoBiTec Ni-IDA is an optimized matrix with low density of IDA ligands. This non-saturating surface concentration of IDA eliminates almost all non-specific interactions of contaminating host proteins with the adsorbent. As a result, MoBiTec Ni-IDA provides higher target protein purity.

IDA is a tridentate chelator which occupies three of the six binding sites in the coordination sphere of the  $\text{Ni}^{2+}$  ion. The remaining three coordination sites are usually occupied by water molecules and can be exchanged with histidine residues of the recombinant protein.

Reliably test your experimental protocols by using MoBiTec's His<sub>6</sub>-GFP bacterial expression vector. His<sub>6</sub>-GFP protein - extracted from *E. coli* - can be detected by eye when loaded on, washed on, and eluted from the Ni-IDA columns. Every step can be quantified easily.

**Ask for  
Free Samples!**

### Order Information, Shipping and Storage

Order No.	Description	Amount
PR-HTK105	MobiSpin Ni-IDA Columns	5 Columns
PR-HTK110	MobiSpin Ni-IDA Columns	10 Columns
PR-HTK004	MoBiTec Ni-IDA Columns	4 Columns
PR-HTK010	MoBiTec Ni-IDA Columns	10 Columns
TOPO-HIS01	His <sub>6</sub> -GFP Bacterial Expression Vector	10 µg DNA
PEG01*	pEG-His1 Vector	5 µg DNA

Columns contain dry matrix; shipped at RT; storable at RT for at least 1 year.

\*Convenient cloning and subsequent purification of your protein of interest: pEG-His1 vector for IPTG inducible expression of gene products with C-terminal RGS motif and His<sub>6</sub>-tag. 5' and 3' primers for sequencing are provided in addition.



# Chapter V

## **Compact Reaction Columns (CRC) Mobicols with Enzymatically Active Matrices**

CRC with Immobilized Proteases

CRC with Immobilized Nucleases

CRC with other Immobilized Enzymes

## Immobilized Enzymes for solid Phase Enzyme Reactions

Compact Reaction Columns (CRC) are small volume columns (Mobicols) containing a matrix with covalently immobilized enzymes. The patented immobilization chemistry was developed in cooperation with the Max Planck Society. It enables us to offer enzymes bound to F7m and G3m matrices with high enzymatic densities. The high enzyme densities of the CRCs result in fast substrate to product turnover. Enzyme reaction occurs when the substrate is loaded onto the column. The sample is recovered from the column quantitatively either by elution or by centrifugation. Since the enzymes are covalently bound, they remain in the column after reaction and product elution, making the CRC reusable. The eluted product solution is free of enzyme and does not require enzyme removal steps. This results in shorter lab protocols and handling procedures.

### Features

- Covalently bound enzymes
- No product contamination
- Suitable for small and large volumes
- Luer-lock connection to reservoir possible
- Short reaction protocol
- CRC is reusable
- Simple product recovery by centrifugation or elution
- Convenient washing enabled/possible
- High enzymatic densities

### Product Description


The immobilized enzymes are stable in aqueous media at a pH range of 5 to 10 and column bleeding is negligible. The “stiff” linkers, which keep the enzymes from the matrix surface, effectively eliminate steric hindrance. This results in high activity of enzymes in the immobilized state.

Column elution characteristics vary depending on the nature of the matrix used for immobilization.

## Cut down time and costs

with immobilized enzymes

MoBiTec

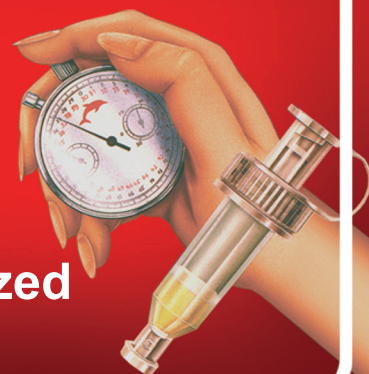

- Columns with enzymes immobilized on matrix F7m are designed for applications involving either large substrate quantities, large volumes and/or large molecules.
- Columns with enzymes immobilized on G3m matrix are designed for applications involving small substrate quantities, small volumes and/or small molecules. Samples with a molecular weight  $\geq 103$  Daltons are recovered in the void volume when loaded on these columns.

This well established technology allows you to expose your reaction solutions to very high concentrations of modifying enzymes. The exposure to high enzyme concentrations allows very short reaction times.

### Compact Reaction Columns (CRC)

Immobilized enzymes are supplied in versatile compact reaction columns (CRC) which fit into 1.5 or 2 ml microcentrifuge tubes. The columns have Luer-lock fittings, allowing direct syringe application of substrate solution, continuous flow processing of bulk solutions, or application of pressure for recovering the substrate. For small substrate volumes (approx. 50  $\mu$ l or less), most enzyme columns can be spun dry in benchtop centrifuges for fast, effective recovery.






## Procedure

- Equilibrate column with reaction buffer by flow-through (since column is provided with storage buffer).
- Load substrate solution in reaction buffer.
- Incubate.
- Elute product solution by:
  - a) centrifugation (smaller volumes, matrix G3m)
  - b) continuous flow mode (larger volumes, matrix F7m)

## Advantages


- Very high enzyme activities
- No enzyme removal
- G3m: no sample dilution
- No long lasting gravity flow
- No chromatography system
- Easy handling, multiple sample processing in a few minutes

**Easy handling,  
brief procedure  
in 4 steps**

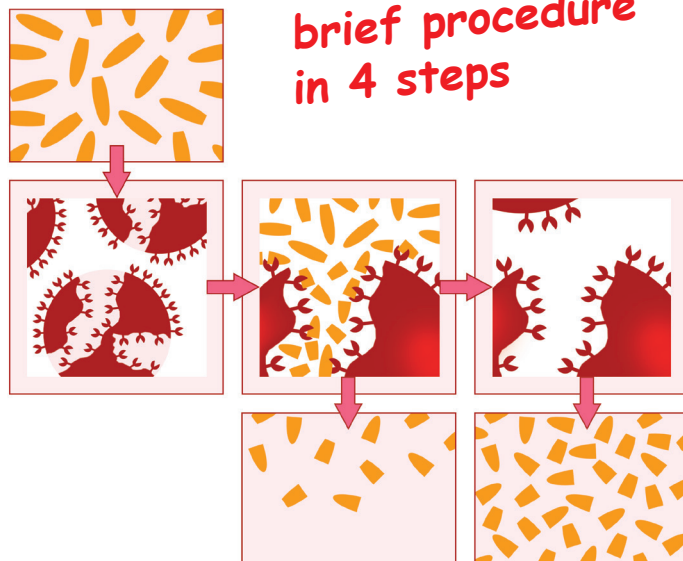

## Matrix Characteristics

Two different immobilization matrices are available: F7m and G3m. They have different pore characteristics, as described below.

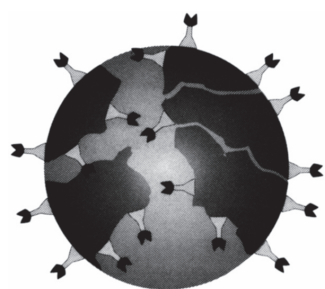


**F7m**  
**Large pores-**  
**high activity!**

### Matrix F7m

Matrix F7m has large pores. Molecules with up to  $10^7$  Dalton molecular weight (most enzymes and substrates) can enter these pores. The total surface of the material (including the surface inside the pores) is very large, resulting in an extremely high enzyme activity on the matrix. The product, which also enters the pores, has to be washed out of the matrix, leading to a dilution of the sample solution by additional 200  $\mu$ l to 500  $\mu$ l, depending on the substrate. If a dilution of the sample solution by 200  $\mu$ l or more cannot be tolerated, matrix G3m should be used.

### Matrix G3m

G3m has small pores; this results in excellent recovery characteristics (i.e. complete recovery in very small volumes). Molecules larger than  $10^3$  Dalton molecular weight (larger peptides, proteins and nucleic acids) cannot enter these pores. The total surface area of the material (only the outside of the spheres) is smaller than for F7m resulting in a smaller enzyme activity on the matrix. The immobilized enzymes as well as most substrates reside only in the volume between the matrix beads (in the dead volume). Thus, the sample solution can be completely removed by centrifugation without dilution. The design of the CR columns makes centrifugation handling easy. Sample solutions down to 20  $\mu$ l can be treated and recovered.

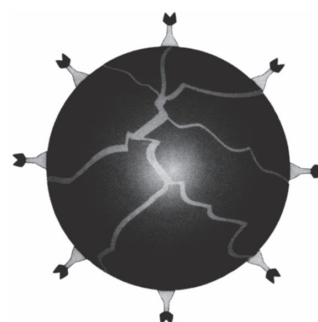

**G3m**  
**Small pores-**  
**no sample dilution!**




## Immobilization Chemistry

The immobilization chemistry has the particular advantage of being extremely stable in aqueous solution in the pH range 5 to 10. Bleeding-out is negligible. The rather stiff linkers separate the enzymes from the matrix surface thus eliminating steric hindrance. In the immobilized state the enzymes retain a high activity.

### Bleeding out in Comparison to Other Immobilization Chemistries

1,000 units alkaline phosphatase were immobilized on the same volume of different activated resins. 200  $\mu$ l Tris/HCl buffer pH 8.0 were applied on the columns for 20 minutes and recovered. The activity of the alkaline phosphatase in the sample was assayed. The relative amounts are shown in the figure below.

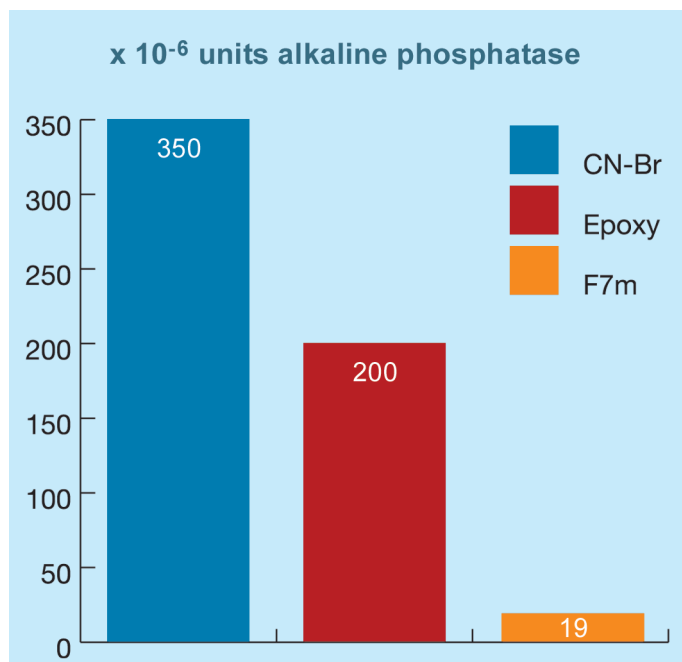

Comparison of three different matrices. According to these results, F7m is the matrix of choice, since bleeding out is significantly reduced compared to CNBr or Epoxy matrices.

## Columns, Filters and Buffers

Our 1 ml Mobicol column is specifically designed for this purpose. The enzymes are immobilized on 200  $\mu$ l matrix which is placed between two filters. The lower filters in the G3m columns have a pore size of 35  $\mu$ m, in F7m columns the lower filters have 10  $\mu$ m pore size. All upper filters have 90  $\mu$ m pore size. The filter material is of polyethylene. The enzyme matrix is delivered in storage buffer. Reaction, washing and storage buffers are provided for the first applications. All buffers are specified on the detailed data sheets provided with the columns.

The amount of liquid, which can be spun out of the CR-column to be replaced by new solution, is 90  $\mu$ l for F7m matrix and 76  $\mu$ l for G3m matrix.

### Effect of pH

The protein coupling to the matrix is stable in the range of pH 5 to pH 10. For short time periods (up to a few hours) the columns can be used also in the range pH 4 to pH 11. The immobilized enzymes should not, however, be stored below pH 5 or above pH 10. At extreme pH values the enzyme/matrix binding will degenerate.

### Biological Activity

DNA plasmids are biologically as active (in transformation) after passage through a CR-column (Proteinase K, RNase A) as the untreated molecules.






## Applications

### Removal of genomic DNA


Immobilized DNase I on G3m matrix is a useful tool for the removal of contaminating genomic DNA during RNA purification for single-cell cDNA-PCR: Standard DNase I digestion of RNA prepared from a single cell often leads to loss of material due to the required extraction step. Therefore, a method eliminating phenol extraction is extremely helpful. Using immobilized DNase I instead of soluble enzyme allows the efficient amplification of cDNA by PCR at the single-cell level in the absence of contaminating DNA. PCR is a process covered by patents owned by Hoffmann La-Roche.

### DNA digestion in a column

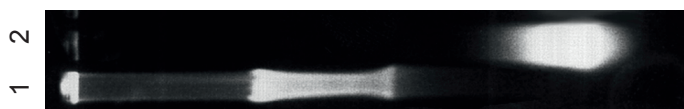

Immobilized DNase I on matrix G3m Agarose gel showing DNA before (lane 1) and after (lane 2) treatment on a DNase I column. The DNA is completely digested on the column.

### Easy digesting RNA during plasmid preps without phenol extractions

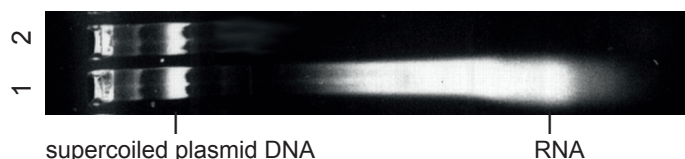


Agarose gel after separation of a clear lysate solution of a plasmid preparation containing DNA (supercoiled and nicked) and RNA (lane 1). Lane 2: Same solution after passage through a RNase G3m column. While the DNA remains unchanged the RNA is digested and cannot be detected.

### Linear vector DNA dephosphorylation using an alkaline phosphatase G3m CR-column

When cloning DNA fragments into plasmids such as expression vectors (i.e. plasmids up to 10 kb), the background resulting from recircularized plasmid

vectors has to be minimized. This can be done by an alkaline phosphatase treatment of the linearized plasmid. Both ends will be dephosphorylated, thus preventing recircularization. The DNA fragments to be cloned, generated by restriction enzyme digestion, will provide the 5' terminal phosphate and thus ligate efficiently to the vector. Following plasmid dephosphorylation, the phosphatase activity must be removed before ligation in order to prevent dephosphorylation of the insert. To overcome this problem CR-columns with immobilized alkaline phosphatase provide dephosphorylated DNA ready for ligation without phosphatase activity contamination.

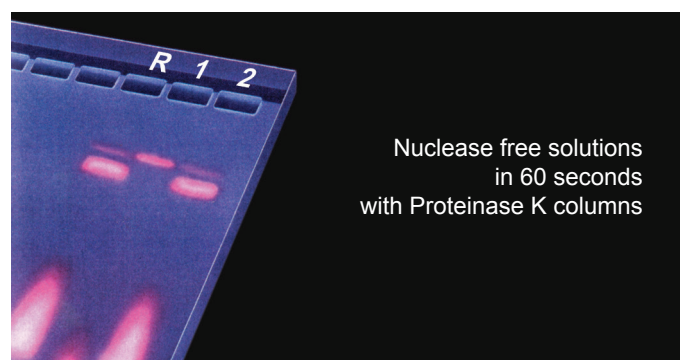

A DNA and RNA solution (reference, lane R) is contaminated with nucleases. After passage through a Proteinase K CR column, both solutions are incubated at 37 °C for one hour, followed by agarose gel analysis. Untreated solution is digested (lane 1) while the recovered sample from the column (lane 2) is free of degradation.

### Insulin Box digestion by an endoprotease Glu-C Compact Reaction Column (matrix G3m)

Since the endoprotease is immobilized, it is completely retained in the column. Thus, the cleaved protein is entirely free of proteinases (and protease fragments); there is no need for further purification of the product.



## Order Information, Shipping & Storage

Order No.	Description	Amount
<b>Immobilized Proteinases (in CRC)</b>		
<b>P5101</b>	Endoproteinase Glu-C (Prot.V8) (one column with 200 µl matrix F7m)	900 U
<b>P3102</b>	Endoproteinase Glu-C (Prot.V8) (one column with 200µl matrix G3m)	22 U
<b>P5401</b>	Papain (one column with 200 µl matrix F7m)	23 U
<b>P3402</b>	Papain (one column with 200 µl matrix G3m)	0.6 U
<b>P5121</b>	Pepsin (one column with 200 µl matrix F7m)	16 mAnson U
<b>P3122</b>	Pepsin (one column with 200 µl matrix G3m)	0.4 mAnson U
<b>P5501</b>	Proteinase K (one column with 200 µl matrix F7m)	27 mAnson U
<b>P3502</b>	Proteinase K (one column with 200 µl matrix G3m)	0.7 mAnson U
<b>P5301</b>	TLCK-α-Chymotrypsin (one column with 200 µl matrix F7m)	55 U
<b>P3302</b>	TLCK-α-Chymotrypsin (one column with 200 µl matrix G3m)	1.4 U
<b>P5701</b>	TPCK-Trypsin (one column with 200 µl matrix F7m)	10,200 St-U
<b>P3702</b>	TPCK-Trypsin (one column with 200 µl matrix G3m)	260 St-U

### Immobilized Nucleases in CRC (except N3403)

<b>N5401</b>	DNase I (one column with 200 µl matrix F7m)	3,500 U
<b>N3402</b>	DNase I (one column with 200 µl matrix G3m)	88 U
<b>N3403</b>	DNase I (kit with 200 µl matrix and 5 empty columns)	88 U
<b>N5101</b>	RNase A (one column with 200 µl matrix F7m)	50 Kunitz-U
<b>N3102</b>	RNase A (one column with 200 µl matrix G3m)	2.5 Kunitz-U

### Immobilized Other Enzymes (in CRC)

<b>A5201</b>	Alkaline phosphatase (CIP) (one column with 200 µl matrix F7m)	1,000 U
<b>A3202</b>	Alkaline phosphatase (CIP) (one column with 200 µl matrix G3m)	100 U
<b>A5101</b>	β-Galactosidase (one column with 200 µl matrix F7m)	600 U
<b>A3102</b>	β-Galactosidase (one column with 200 µl matrix G3m)	15 U

Shipped at RT; store at 4 °C

**Note: Never freeze a CR-column!**

### All Columns are delivered with

- 200 µl immobilized enzyme matrix
- concentrated buffers for the first applications
- our CRC handbook about immobilized enzymes
- a data sheet with the description of the immobilized enzyme, the buffer composition and the protocol for the use of the column.

**A detailed handbook including all protocols is provided with the product. It is also available for download on our website**

**[www.mobitec.com](http://www.mobitec.com)**



## References

- William N. Addison, Fereshteh Azari, Esben S. Sørensen, Mari T. Kaartinen, and Marc D. McKee Pyrophosphate Inhibits Mineralization of Osteoblast Cultures by Binding to Mineral, Up-regulating Osteopontin, and Inhibiting Alkaline Phosphatase Activity *J. Biol. Chem.*, May 2007; 282: 15872 - 15883.
- Jae Min Chong and David W. Speicher Determination of Disulfide Bond Assignments and N-Glycosylation Sites of the Human Gastrointestinal Carcinoma Antigen GA733-2 (CO17-1A, EGP, KS1-4, KSA, and Ep-CAM) *J. Biol. Chem.*, Feb 2001; 276: 5804 - 5813.
- Scott Ficarro, Olga Chertihin, V. Anne Westbrook, Forest White, Friederike Jayes, Petr Kalab, Jarrod A. Marto, Jeffrey Shabanowitz, John C. Herr, Donald F. Hunt, and Pablo E. Visconti Phosphoproteome Analysis of Capacitated Human Sperm. Evidence of Tyrosine Phosphorylation of a Kinase-Anchoring Protein 3 and Valosin-Containing Protein/P97 During Capacitation *J. Biol. Chem.*, Mar 2003; 278: 11579 - 11589.
- Marcin Gierdalski and Sharon L. Juliano Factors Affecting the Morphology of Radial Glia, Cereb Cortex, Jun 2003; 13: 572 - 579.
- Jeff Holderness, Larissa Jackiw, Emily Kimmel, Hannah Kerns, Miranda Radke, Jodi F. Hedges, Charles Petrie, Patrick McCurley, Pati M. Glee, Aiyappa Palecanda, and Mark A. Jutila Select Plant Tannins Induce IL-2R $\alpha$  Up-Regulation and Augment Cell Division in  $\gamma\delta$  T Cells *J. Immunol.*, Nov 2007; 179: 6468 - 6478.
- R.E. Leube, U. Leimer, C. Grund, W.W. Franke, N. Harth, and B. Wiedenmann Sorting of synaptophysin into special vesicles in nonneuroendocrine epithelial cells *J. Cell Biol.*, Dec 1994; 127: 1589.
- Thomas S. Nühse, Thomas Boller, and Scott C. Peck A Plasma Membrane Syntaxin Is Phosphorylated in Response to the Bacterial Elicitor Flagellin *J. Biol. Chem.*, Nov 2003; 278: 45248 - 45254.
- Thomas S. Nühse, Allan Stensballe, Ole N. Jensen, and Scott C. Peck Large-scale Analysis of in Vivo Phosphorylated Membrane Proteins by Immobilized Metal Ion Affinity Chromatography and Mass Spectrometry *Mol. Cell. Proteomics*, Nov 2003; 2: 1234.
- A. Pries, H. Priefert, N. Krüger, and A. Steinbüchel. Identification and characterization of two *Alcaligenes eutrophus* gene loci relevant to the poly( $\beta$ -hydroxybutyric acid)-leaky phenotype which exhibit homology to ptsH and ptsI of *Escherichia coli*. *J. Bacteriol.*, Sep 1991; 173: 5843 - 5853.
- J. Ronnenberg, B. Preitz, G. Wostemeier, and S. Diekmann Immobilized residue-specific endoproteases for protein sequencing. *J Chromatogr B Biomed Appl*, Jun 1994; 656(1): 169-77.
- N. M. Sawtell, D. K. Poon, C. S. Tansky, and R. L. Thompson The Latent Herpes Simplex Virus Type 1 Genome Copy Number in Individual Neurons Is Virus Strain Specific and Correlates with Reactivation *J. Virol.*, Jul 1998; 72: 5343 - 5350.
- Richard L. Thompson and N. M. Sawtell Pathogenesis and Immunity: Replication of Herpes Simplex Virus Type 1 within Trigeminal Ganglia Is Required for High Frequency but Not High Viral Genome Copy Number Latency *J. Virol.*, Jan 2000; 74: 965 - 974.
- R. L. Thompson and N. M. Sawtell Pathogenesis and Immunity: Herpes Simplex Virus Type 1 Latency-Associated Transcript Gene Promotes Neuronal Survival *J. Virol.*, Jul 2001; 75: 6660 - 6675.





# Fax Back - Order Form

**Please complete this form and  
fax back to: + 49 551 70722-22**

Your local distributor:

## Delivery Address

First name \_\_\_\_\_

Last name \_\_\_\_\_

Institution/Company \_\_\_\_\_

Street, No. \_\_\_\_\_

ZIP, City \_\_\_\_\_

Country \_\_\_\_\_

Phone \_\_\_\_\_

Fax \_\_\_\_\_

Email \_\_\_\_\_

## Billing Address (if different from delivery address)

First name \_\_\_\_\_

Last name \_\_\_\_\_

Institution/Company \_\_\_\_\_

Street, No. \_\_\_\_\_

ZIP, City \_\_\_\_\_

Country \_\_\_\_\_

Phone \_\_\_\_\_

Fax \_\_\_\_\_


Email \_\_\_\_\_

Order #	Description	Quantity	Unit	Price

Date \_\_\_\_\_

Signature \_\_\_\_\_



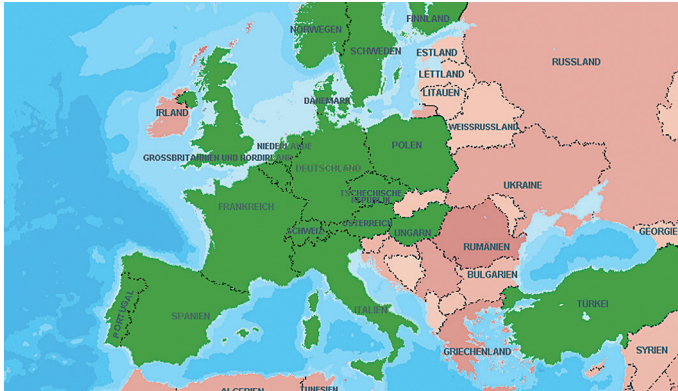

Do you seek a  
**powerful partner**  
to market your ideas  
professionally worldwide?




#### Approach MoBiTec today with

- ideas that are ready for commercialization,
- fully developed products and technologies for licensing,
- technology concepts for further development,
- proposals for scientific collaborations.

New technology offers will be assessed, from both a technical and commercial perspective, by our R&D and commercial groups.

Therefore, we need appropriate information for an initial review of the technology. Further details of the technology may be required to reach a final conclusion. The assessment of your technology will be performed under the protection of a confidentiality agreement.

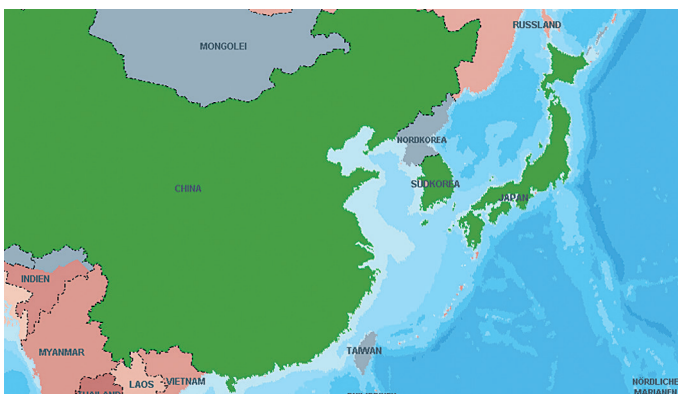

#### MoBiTec GmbH

Phone +49 551-70722-0  
Fax +49 551-70722-22  
E-mail [info@mobitec.de](mailto:info@mobitec.de)  
Internet [www.mobitec.com](http://www.mobitec.com)

Visit us at <http://www.mobitec.com>





# Mobicols

SEPARATION

SEPARATION

ISOLATION

ISOLATION

TREATMENT

TREATMENT

PURIFICATION

PURIFICATION

Your local distributor:

[www.mobitec.com](http://www.mobitec.com)

MoBiTec GmbH,  
Lotzestrasse 22a, 37083 Goettingen, Germany

Phone: +49 (0)551-70722-0

Fax: +49 (0)551-70722-22

E-mail: [info@mobitec.com](mailto:info@mobitec.com)



**Mo Bi Tec**  
MOLECULAR BIOTECHNOLOGY