

Xpert Fast Hotstart DNA Polymerase (dNTPs included)

#GE25.0250 (250 U) | GE25.1000 (1000 U) | GE25s (trial size)
(FOR RESEARCH ONLY)



Product: Xpert Fast Hotstart DNA polymerase is a robust enzyme, suitable for the amplification of difficult targets with extreme speed (4-8kb/min), yield and specificity. The buffer composition has been optimized via high-throughput screening and makes the enzyme particularly resistant to PCR inhibitors. As a result, Xpert Fast Hotstart DNA Polymerase is suitable for direct PCR from unpurified samples, including blood, urine, and bacterial colonies, making this the ideal choice for consistent results in fast complex PCR amplifications

Applications: Direct PCR, Fast amplification of complex targets (including GC-rich and AT-rich templates), Colony PCR, Multiplex PCR.

Contents: Xpert Fast Hotstart DNA Polymerase is supplied at a concentration of 5U/μl. The reaction buffer (5x) includes not only 15 mM MgCl₂, but also 5mM dNTPs, enhancers and stabilizers, optimized to increased PCR success rates.

Component	GE25.0250	GE25.1000	GE25s
Xpert Fast Hotstart (5U/μl)	50 μl	4x 50 μl	50 μl
Xpert Fast Hotstart Reaction Buffer (5x), containing dNTPs	2x 1ml	8x 1ml	400 μl

QC: Functionally tested in PCR. Absence of endonucleases, exonucleases, and ribonucleases was confirmed by appropriate assays

Properties:

- Amplicon size: up to 5kb
- Extension rate: 2sec/kb (for targets up to 1kb) and 15sec/kb for larger targets
- Hotstart: Yes
- A-overhang: Yes

Storage: -20°C and protected from light for at least 1 year. No loss of performance is detected after 20 freeze/thaw cycles.

Prior to use

Optimal PCR cycling conditions (incubation times, temperatures, and concentrations) depend on DNA target (GC-content, size, quantity, purity, etc), specific primers, buffer composition, MgCl₂-concentration, and Xpert Fast Hotstart-concentration and need to be determined case by case. Xpert Fast Hotstart Reaction Buffer includes dNTPs and has been already been optimized with respect to the MgCl₂-concentration and other components to maximize success rates. It is not recommended to add additional MgCl₂ or other PCR enhancers. We suggest to start with the following basic protocol and subsequently optimize annealing temperature, incubation times and cycling number.

Basic Protocol

1. Mix for each PCR reaction, starting with the greatest volume (usually water) and ending with Xpert Fast:

Component	Volume (25µl)	Final Conc.
Xpert Fast Hotstart Reaction Buffer (5x)	5 µl	1X
MgCl ₂ (already included in PCR Buffer)	0 µl	3 mM
dNTPs (already included in PCR Buffer)	0 µl	1 mM (each)
Forward primer (5 pmol/µl)	2 µl	0.4 µM
Reverse primer (5 pmol/µl)	2 µl	0.4 µM
Template DNA*	0.25 – 10 µl	1-250 ng*
Xpert Fast Hotstart DNA Polymerase (5U/µl)	0.1 µl**	0.5 U**
PCR-grade water	up to 25 µl	

*) In case of cDNA <50ng and in case of gDNA <250ng (total amount).


In case of Direct PCR: use 1µl of mammalian blood or urine per 25µl reaction volume

In case of Colony PCR: use a sterile tip to resuspend a colony (or 5µl of an overnight culture) directly in the PCR reaction mix.

***) For very difficult targets, e.g. higher enzyme concentration may be required (up to 2.5U)

In order to minimize risk of contamination, reagent loss and improve pipetting accuracy, we recommend to prepare a mastermix for multiple samples (N), always including a negative control for the detection of possible contaminants, by mixing all components (N+1), except template DNA, dividing the mixture equally into each tube and then add template DNA or PCR grade water in case of the control to the individual PCR tubes.

2. Set-up initial PCR amplification as follows:

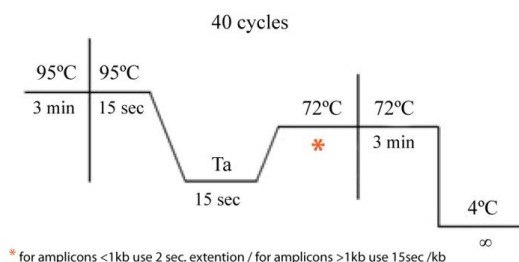

After an initial cycle of 3 min at 95°C (enzyme activation, denaturation of template DNA), cycle 40 times for 15 seconds at 95°C, 15 seconds at Ta, and 2 to 75 seconds at 72°C for extension (2 seconds for targets below 1 kb and 15 seconds per kb for target DNA up to 5kb | or 90 seconds total per cycle in case of multiplex PCR, independent of the sizes of the amplicons). Set the annealing temperature (Ta) as the melting temperature (Tm) of the primer with the lowest Tm. Include a final extension step of 3 min at 72°C to ensure that all amplicons are fully extended and include 3'-A-overhang. Analyze PCR products by DNA Agarose gel electrophoresis.

NOTE: In case of Colony PCR, increase initial denaturation time to 10 minutes.

Optimization

Annealing Temperature (Ta) and Primers

Optimizing the annealing temperature is crucial, especially in case of multiplex PCR, as a too low temperature might result in non-specific amplification whereas a too high temperature results in no amplification. The melting temperature (Tm) is defined as the temperature in which 50% of the primer and its complementary sequence of the target DNA are present as duplex DNA. By increasing the temperature above the melting temperature, this percentage decreases, however, primers will still anneal (up to a certain point) and initiate extension. PCR can therefore be performed at temperatures of several degrees higher than Tm and it is therefore recommended to optimize the Ta by performing a temperature gradient (e.g. starting at the lowest Tm or a few degrees below and increasing with 2°C increments). Ideally, primers have melting temperatures of approximately 60°C and final concentration should be between 0.2 and 0.6µM (each).

Incubation times and number of cycles.

Denaturation and annealing steps may require less time depending on the thermocycler apparatus (ramp rate), reaction volume and PCR tube (varies with the efficiency of heat-transfer). It might be worthwhile to optimize (reduce) times to as low as 10 seconds for both denaturation and annealing steps, which will greatly reduce overall PCR time. **Fast cycling conditions should not be applied in multiplex PCR.** The initial recommended 90 seconds per cycle may be further extended to increase yield. It might be worthwhile to reduce number of cycles from 40 to 25-30, depending on the success of amplification.