

Cloning Vector p3T

**Order #
P123T**



Mo Bi Tec
MOLECULAR BIOTECHNOLOGY



Contents

1. Summary	3
2. Introduction	3
3. The p3T Vector	4
4. Protocols	5
4.1. Amplification of DNA	5
4.2. Cloning by Polyadenylation	5
5. Buffers	7
6. Bacterial Host	7
7. References	7
8. Order Information, Shipping and Storage	8



1. Summary

The p3T vector provides a flexible system for the direct cloning of PCR products. Based on pBluescript II SK + it allows the cloning of PCR products via both single and multiple T overhangs. Utilizing a unique series of restriction sites the same vector can be cleaved by different enzymes to produce 1, 2 or 3 T overhangs. This permits either the direct cloning of PCR products (via a single A extension) or the option of polyadenylating the PCR fragment and cloning via multiple A extensions. This procedure has the advantage that less material is required for cloning. The vector also contains a SmaI site to reduce vector background. MscI sites (blunt) flank the insert allowing excision with minimal flanking regions. Blue/white selection by α -complementation can still be used with the p3T vector.

2. Introduction

A number of vectors exist for the direct cloning of PCR products (Holten, T.A. et al. 1991; Kovalic, D. et al. 1991; Marchuck, D. et al. 1991; Mead, D.A. et al. 1991). Many of these are based on the addition of single T residues or tandem restriction sites that can be cleaved excising a small fragment to produce single T-overhangs. A major disadvantage is that due to the inefficient ligation of single base overhangs a large amount of PCR fragment is required and/or the use of extremely efficient competent cells. The p3T vector (Mitchell, D.B. et al., 1992) can be cut with PflMI to generate a 3-T overhang. If the PCR fragment is polyadenylated using terminal deoxynucleotidyl transferase available at MoBiTec, see chapter 8, it can be cloned with high efficiency. The polyadenylation is a simple procedure requiring only a five minute reaction time.

The p3T vector was constructed using a pBluescript II SK+ backbone. The cassette was synthesized and cloned between the PstI and BamHI sites in the polylinker. The insertion is 48 nucleotides long and preserves the β -galactosidase reading frame permitting blue/white selection. A naturally occurring BcgI site in the ampicillin gene was removed by site directed mutagenesis. This vector is shown in Figure 1.

PCR products can be cloned directly using XcmI cut vector or, more efficiently, cloned using PflMI cut vector and polyadenylating the PCR fragment. XcmI leaves a single dT overhang, whereas PflMI produces a multiple dT overhang (cleavage sites see chapter 4.2.1.).

Advantages

- single or multiple T-overhangs can be utilized
- direct cloning of PCR products (single dA extension)
- or polyadenylation and cloning via multiple dA extensions
- requires less amplified DNA
- SmaI site to reduce the vector background
- MscI sites flank the insert for optimal excision
- cloning with high efficiency
- blue/white selection by α -complementation



3. The p3T Vector

A schematic presentation of the p3T vector is shown in Figure 1. The complete sequence is deposited in the EMBL database (accession number Z46733) or can be obtained from MoBiTec by E-mail.


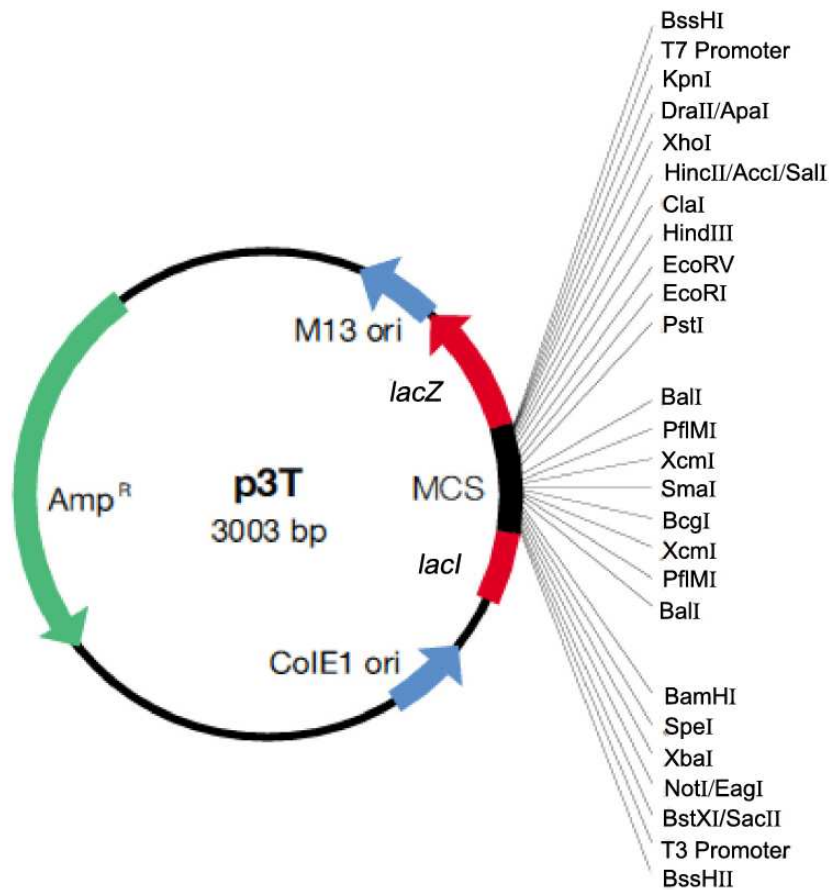


Figure 1. The p3T vector.

Amp^R: ampicillin resistance gene;
ColE1 ori: pBR322 derived origin of replication;
lacI: β -galactosidase promoter sequence;
lacZ: the alpha peptide of β -galactosidase;
fd: the f1 filamentous phage derived origin of replication.

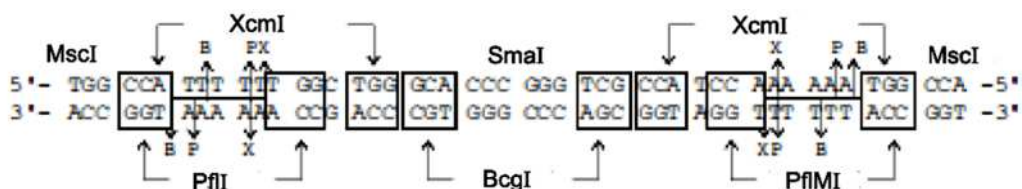

The insert sequence in p3T.

The insert orientation in the plasmid is from PstI to BamHI (left to right). The restriction enzyme recognition sites are boxed and the sites of cleavage are indicated by arrows. B, BcgI; P, PfiMI; X, XcmI.



4. Protocols

4.1. Amplification of DNA

Use standard reaction conditions.

Note: *MoBiTec tubes for thermocyclers do not require oil.*

4.2. Cloning by Polyadenylation

- For a single dT overhang, digest 5 µg of p3T with 20 - 30 U of restriction enzyme *Xcm* I for about 4 hours at 37 °C. The enzyme concentration and incubation time may be reduced.

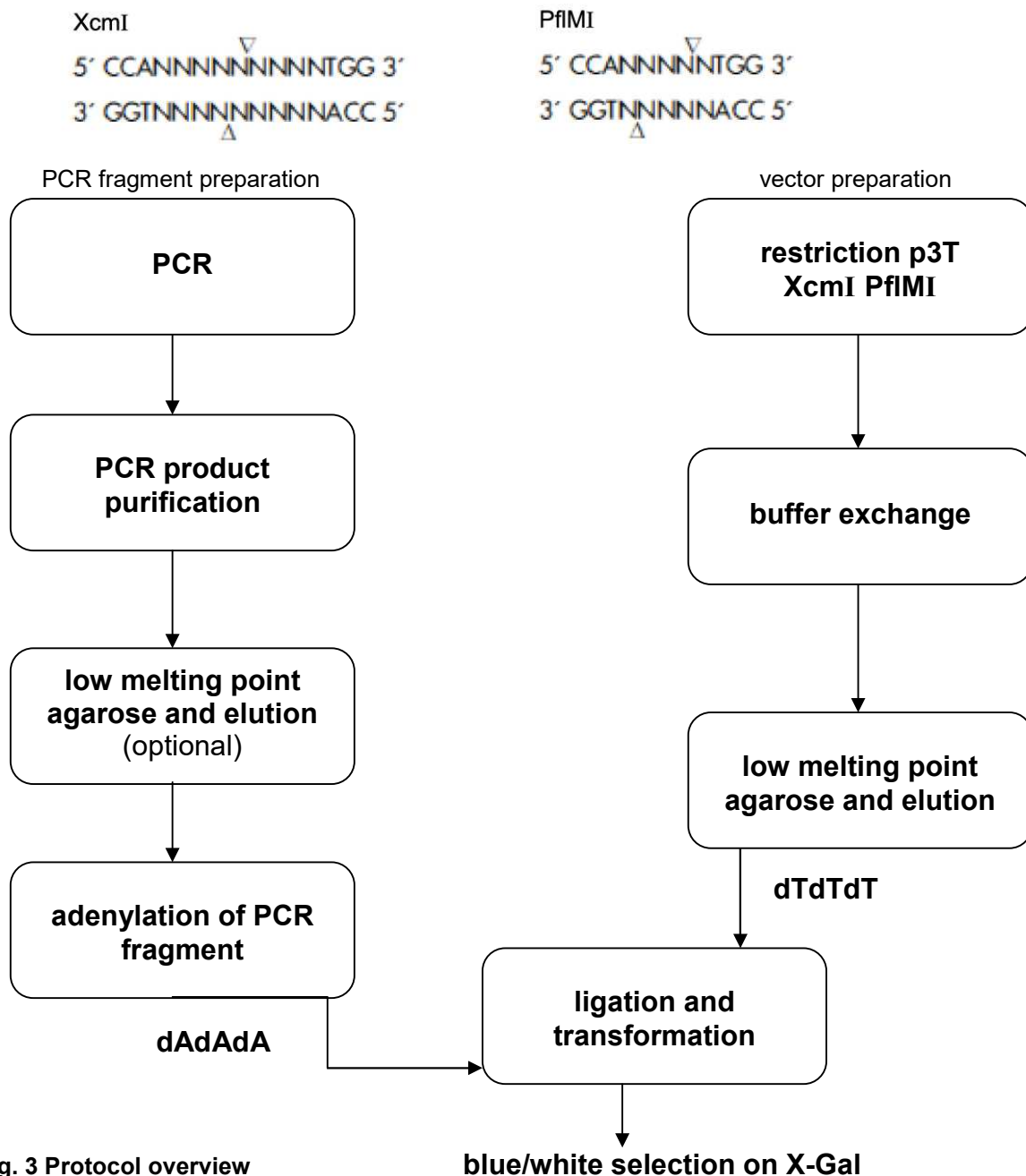


Fig. 3 Protocol overview

blue/white selection on X-Gal



- For a multiple dT overhang digest 5 µg of p3T with 20 - 30 U of restriction enzyme *Pfl* M1 for about 4 hours at 37 °C. The enzyme concentration and incubation time may be reduced.
- Confirm by analysis on a 0.8% agarose gel, that the digestion is complete.
- Precipitate DNA with 1/10 vol 3 M sodium acetate pH 5.2 and 2 volumes of 100% ethanol.
- Resuspend the pellet in an appropriate volume of 1x TE pH 7.5 and loading buffer, run on a 0.8% low melting point agarose gel and elute, in order to separate the small excised fragment from the linearized vector.

Note: *To simplify the elution procedure, the agarose with the appropriate band of DNA may be frozen and then filled into an empty Mobicol with filter (see chapter 8). After mixing the agarose with buffer, the Mobicol is centrifuged within an Eppendorf tube. The agarose remains on the filter and most of the DNA is eluted into the Eppendorf tube. An ethanol precipitation might be advisable after this procedure.*

4.2.2. Elution and adenylation of PCR fragment

- Extract the PCR amplified DNA with an equal volume of chloroform isoamylalcohol (25:1) and precipitate with 1/10 vol 3 M NaOAc pH 5.2 and 2 volumes of 100 % ethanol.
- Resuspend the pellet in an appropriate volume of 1x TE pH 7.5 and loading buffer, run on a 1.2 - 1.5% low melting point agarose gel and recover the fragment. If the produced PCR fragments are very uniform, the agarose gel and elution may be omitted.
- For the ligation to multiple dT overhangs, set up the 50 µl adenylation reaction described below.
- | | |
|----------------|--|
| 10 µl | 5x adenylation buffer (as in chapter 5) |
| 0.1 mM | (final concentration) DTT |
| 1 µl | 1 mM dATP |
| 10 U | terminal deoxynucleotidyltransferase see chapter 8 |
| 0.5 - 1.0 pmol | 3'-OH termini of eluted PCR fragment
(add just before incubation!!) |

H₂O to 50 µl

- Incubate the reaction for 5 minutes. Stop the reaction by adding 2 µl 0.5 M EDTA pH 8.0. Add 50 µl 1x TE pH 7.5 and extract immediately, twice with phenol/chloroform/isoamyl alcohol (25:24:1) and once with chloroform/isoamyl alcohol (25:1). Precipitate with a 1/10 vol 3 M sodium acetate pH 7.0 and 2 volumes of 100% ethanol. Resuspend in 1x TE pH 7.5 at the desired concentration.

4.2.3. Ligation and transformation

Ligation and transformation should be performed according to standard procedures (Sambrook *et al.*, 1989).

A *Sma*I digest after ligation and before transformation decreases the vector background, since it linearizes uncut vectors without insert. If there is the possibility of a *Sma*I site within the insert, this step should be omitted.

Transformants are selected on LB-agar plates containing ampicillin, X-Gal and IPTG (standard protocols).



5. Buffers

TE 10x

100 mM Tris
10 mM EDTA
pH 7.5

5x adenylation buffer as supplied with MoBiTec's Terminal Deoxynucleotidyl Transferase:

500 mM K Cacodylate buffer pH 7.2
10 mM CoCl₂
10 mM DTT (supplied in an extra vial)

6. Bacterial Strains

E. coli standard strains such as XL-1 Blue may be used for propagation of p3T and the selection of transformants.

7. References

Holton, T.A. and Graham, M.W. (1991) *Nucleic Acids Res.* 19 1156.

Kovalic, D., *et al.* (1991) *Nucleic Acids Res.* 19 4560.

Marchuck, D. *et al.* (1991) *Nucleic Acids Res.* 19 1154.

Mead, D.A. *et al.* (1991) *Bio/Technology* 9 657-663.

Mitchell, D.B. *et al.* (1992) *PCR Meth. App.* 2 81-82.

Sambrook *et al.* *Molecular Cloning*, (1989) Cold Spring Harbour

Sandhu, G.S. *et al.* (1989) *Biotechniques* 7 689-690.



8. Order Information, Shipping and Storage

Order#	Product	Quantity
P123T	p3T, PCR Product Cloning Vector, lyophilized DNA	5 µg
shipped at RT; store at 4 °C		
M1002	Mobicols with 2 different screw caps without filters	50
shipped at RT; store at RT		

We Bring The World Of Biotechnology To You



www.bocascientific.com

(781) 686-1631

info@bocascientific.com