

Mobicols with 2 different screw caps without filters


Product Information Sheet
M1002



SUMMARY

shipped at RT; store at RT

For research use only

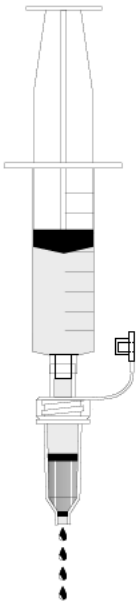

Syringe connected to Luer-lock cap


Mobicol "Classic" is a 1 ml column with particular handling advantages.

Mobicol "Classic"

- is compatible to the laboratory standard,
- can be centrifuged in minifuges,
- has a Luer-lock adapter.
- Smallest and largest volumes can be treated easily.



Filters of the pore sizes 10 μm , 35 μm or 90 μm are available for Mobicols. The filters can be inserted easily with the help of a plastic tool which is delivered with the columns. The upper large filter fits tightly over the whole length of the column body.

Mobicols can be closed tightly with a screw-on cap and a plug at the outlet. The volume is 800 μl . Mobicols can be autoclaved at 120 $^{\circ}\text{C}$ (with filters at 110 $^{\circ}\text{C}$ for 10 minutes).

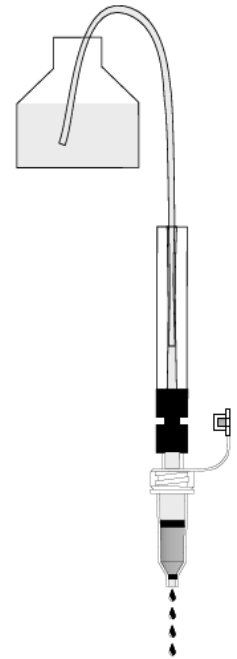

Adapter M3001: tubing to Mobicol

Mobicol "Classic"

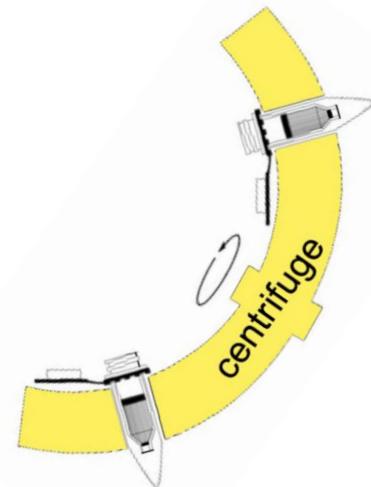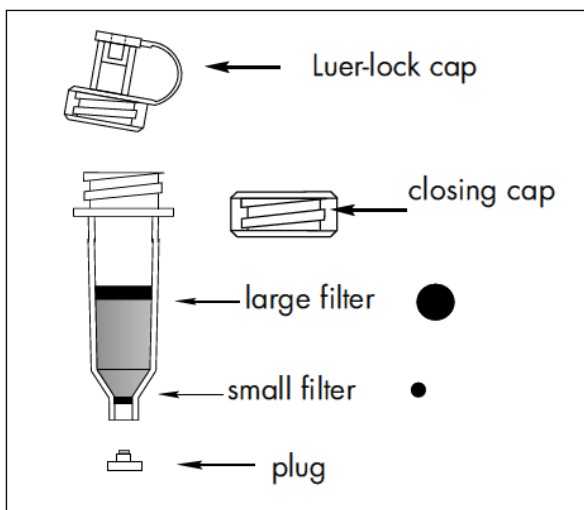

Mobicols can be centrifuged (fit into 1.5 ml or 2.2 ml tubes)

Mobicols with 2 different screw caps without filters

Product Information Sheet # M1002



The Mobicol "Classic" M1002 is supplied with two screw caps (Luer-lock cap and screw cap) and a plug for the outlet. Small and large filters with different diameters (10 μm , 35 μm and 90 μm) are separately available, so you can use the filter of your choice. With this features the Mobicol "Classic" can be used for a variety of applications.

Order Information, Shipping and Storage

Order#	Product	Quantity
M1002	Mobicol "Classic" with 2 different screw caps, without filters	50
M1003	Mobicol "Classic" with screw cap, without filters	50
M2110	Filter (small) 10 μm pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M513515	Filter (small) 35 μm pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M2190	Filter (small) 90 μm pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M2210	Filter (large) 10 μm pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M523515	Filter (large) 35 μm pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M2290	Filter (large) 90 μm pore size, for Mobicol M1002, M1003, M1050 & M1053	50
shipped at RT; store at RT		

Column Accessories

Order#	Product	Quantity
M3001	Tubing-to-Luer-lock adapter (material outside: metal; inside: glass and PTFE; with 1 m PTFE tubing) for Mobicols or laboratory columns	1
M3002	Luer adapter, (M+F) 10 each	20
M3003	Luer adapter for the top-cap (M)	20
M3004	Luer adapter for outlet (F)	20
M3005	Bottom plugs for M1002	50
M3006	Stopcocks, one-way, Luer-lock	4
M3007	Luer-lock adapter (M/M)	1
shipped at RT; store at RT		