

# Mobicol "F" with fixed outlet plug and 1 closed screw cap, without filters



## Product Information Sheet # M1053



### SUMMARY

shipped at RT; store at RT

For research use only

Mobicol "F" is a 1 ml column with particular handling advantages.




Mobicol "F"

- is compatible to the laboratory standard,
- can be centrifuged in microfuges,
- smallest and largest volumes can be treated easily.

Filters of the pore sizes 10  $\mu\text{m}$ , 35  $\mu\text{m}$  or 90  $\mu\text{m}$  are available for Mobicols. The filters can be inserted easily with the help of a plastic tool which is delivered with the columns. The upper large filter fits tightly over the whole length of the column body.

Mobicols can be closed tightly with a screw-on cap and a plug at the outlet. The volume is 800  $\mu\text{l}$  with a small filter and 700  $\mu\text{l}$  with a large filter inserted. Mobicols can be autoclaved at 120  $^{\circ}\text{C}$  (with inserted Filter at 110  $^{\circ}\text{C}$  for 10 min).

### Mobicol "F"

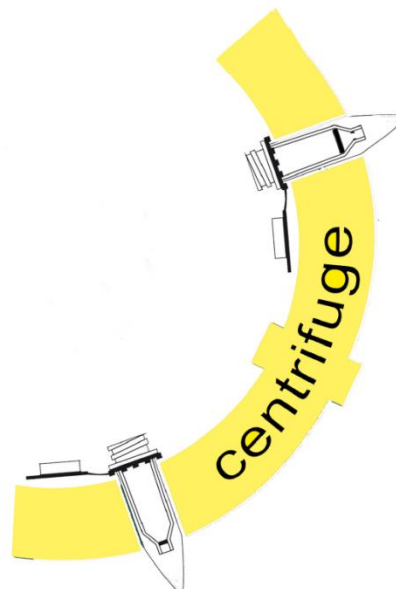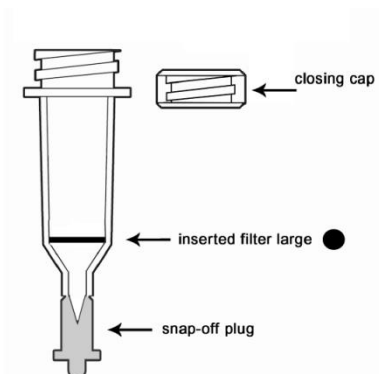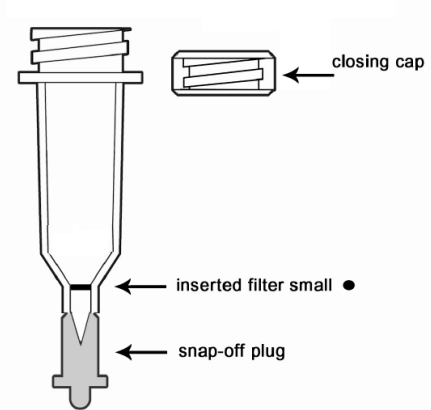


Mobicols can be centrifuged (fit into 1.5 ml or 2 ml tubes)

# Mobicol "F" with fixed outlet plug and 1 closed screw cap, without filters



## Product Information Sheet # M1053

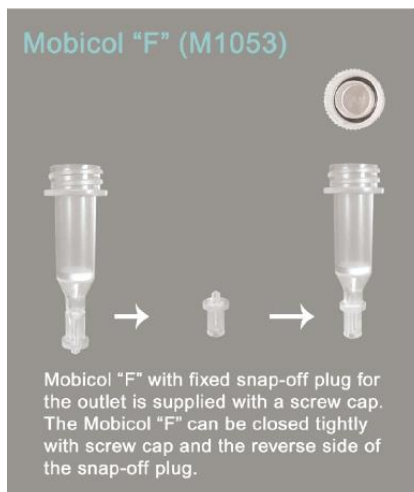

The Mobicol "F" M1053 is supplied with a screw cap and two tools for insertion of filters. The reverse side of the snap-off plug can be used for resealing. Small and large filters with different diameters (10  $\mu\text{m}$ , 35  $\mu\text{m}$  and 90  $\mu\text{m}$ ) are separately available, so you can use the filter of your choice. With this features the Mobicol "F" can be used for a variety of applications.

### Order Information, Shipping and Storage

Order#	Product	Quantity
M1053	Mobicol "F" with fixed outlet plug and screw cap, without filters	50
M2110	Filter (small) 10 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M513515	Filter (small) 35 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M2190	Filter (small) 90 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M2210	Filter (large) 10 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M523515	Filter (large) 35 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
M2290	Filter (large) 90 $\mu\text{m}$ pore size, for Mobicol M1002, M1003, M1050 & M1053	50
shipped at RT; store at RT		